




## TABLES FOR ECHOCARDIOGRAPHY FOR EXAMINATION AND TREATMENT

2021.1



## TABLES FOR SURGERY AND HOSPITAL

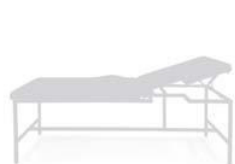

...WITH FIXED HEIGHT



... WITH HEIGHT ADJUSTMENT



# Tables for echocardiography



**Durable** ✓



**Robust** ✓



**Easy-care** ✓




**Cut-out for echocardiography and self fixing upholstery part** ✓




**High load capacity** ✓

Tables with height adjustment additionally:

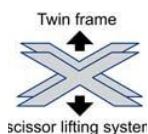

**Twin frame scissor lifting system** ✓

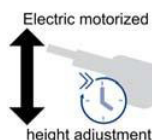

**High adjustment speed** ✓



**Proven Dewert-quality** ✓



## Content

Good reasons for a table for echocardiography ..... 3

Tables for echocardiography with fixed height

2-sectioned tables for echocardiography:

- Model 117 / 177 ..... 6

Tables for echocardiography with height adjustment

Upholstery part(s) positive adjustable

2-sectioned tables for echocardiography:

- Model 2550XLE / 2551XLE ..... 8
- Model 2550XL/H / 2551XL/H ..... 10

3-sectioned tables for echocardiography:

- NEW** • Model 2552XLE ..... 12
- Model 2552XL/H ..... 14

Tables for echocardiography with height adjustment

Upholstery part(s) positive and negative adjustable

- Model 2690XLE ..... 16
- Model 2695XLE ..... 18

Tables for echocardiography with height adjustment and tilt function

- Model 2095XLE ..... 20



All Models with **cut-out for echocardiography**:

- Cut-out for echocardiography and self fixing upholstery part to close
- Conical shape prevents the unintended slipping out of upholstery part




# Tables for echocardiography

## Good reasons for a table for echocardiography




### Models with fixed height

- Frame made of square tube, 30 x 30 x 2 mm
- With cross brace between the feet, longitudinal brace and central support for better stability
- Load capacity: 250 kg

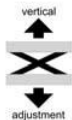



### Electric motorized height adjustment

- Electric motorized height adjustment by AC-motor with a pressure capacity of 10000 N and constant adjustment speed even in case of high load
- EPR-safety power cable with strain relief, length: 4 m
- Integrated safety devices (according to the newest advice of BfArM)
  - Automatic blocking of electromotive drive
  - Active freewheeling downwards, automatic declutch in case of resistance
- Plug connector for potential equalization in acc. to DIN 42801
- Load capacity: 250 kg

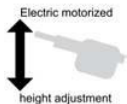

### Hydraulic height adjustment

- Manual foot operated height adjustment by a hydraulic pump with a pressure capacity of 10000 N
- Pedals for foot-controlled height positioning on both sides of the table
- Load capacity: 250 kg

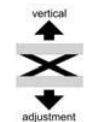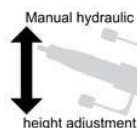

### Dewert-hygiene-equipment

#### Enables optimal cleaning and disinfection

- Artificial leather underneath the upholstery
- Vent holes underneath the upholstery, closed by special caps with valve function
- Easy hygiene because of an open construction: good accessibility of components, no plastic coverings which are breakable in the ageing process
- Option: Innovative artificial leather Skai® Toronto EN with staynu





# Tables for echocardiography

## Good reasons for a table for echocardiography ...



### Upholstery:



- Upholstery thickness: 60 mm, Model 117 / 177: 80 mm
- Multiple glued, high-strength laminated wood panel (thickness: 18 mm) no breaking chipboards in case of pressure on the edges
- Two-layer homogenous polyurethane foam, density 46 and 75, Ökotex-certified
- 2 different variants for artificial leather eligible for upholstery covering:
  - Skai® Pandoria Plus with velvet-smooth surface
  - Skai® Toronto EN with innovative staynu-technology

### Variants:

#### ... with fixed height and 2-sectioned lying surface:



- Two stable parallel mounted safety metal rachets for positive adjustment of head part
- Longitudinal brace for additional stability
- Galvanised steel hinges connect the upholstery parts
- Frame welded as one unit (powder coated) or composed of two units (chromium-plated)

#### ... with height adjustment:



### Technology:



- Extremely stable **twin frame scissor lifting system**, with high torsional stiffness due to thick-walled square tubes with big dimensions, also for under and upper frame and scissor


- High stability due to massive 20 mm axes



- Widely-spaced feet and a high deadweight guarantee stability
- Robust frame construction which is proven in demanding everyday life of hospitals for years



- Abrasion-resistant guide system at under and upper frame due to the use of roller bearings

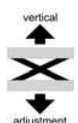



# Tables for echocardiography

## Good reasons for a table for echocardiography ...

### ... with height adjustment and 2-sectioned lying surface:



- Head part with cut-out for echocardiography positive adjustable to 25° by two stable parallel mounted safety metal rachets with tubular steel frame for high torsion rigidity

### ... with height adjustment and 3-sectioned lying surface:

- Head part with cut-out for echocardiography positive adjustable to 25° by two stable parallel mounted safety metal rachets with tubular steel frame for high torsion rigidity
- Fixed middle part
- Foot part positive adjustable by gas spring
- Pressure stable gas spring ensures working at foot part **without vibrations**
- Operating lever to manipulate the gas spring easily accessible



### ... with height adjustment and 3 adjustable upholstery sections:

- Head part stepless positive and negative adjustable by gas spring
- Middle part stepless positive and negative adjustable by electromotive drive  
→ Electromotive drive control by hand switch
- With one or two cut-out(s) for echocardiography and upholstery part(s) to close, depending on model
- Foot part stepless positive and negative adjustable:  
Possible positions including:  
Roof position, Trendelenburg, Reverse Trendelenburg



### ... with height adjustment and tilt function:

- Head part stepless positive adjustable by gas spring
- Independent electric motorized height adjustment and tilt function  
→ Electromotive drive control of tilt function by hand switch

without  
picture

# Tables for echocardiography

## 2-sectioned table for echocardiography with fixed height

... for examination, treatment ...

**Model 117-00 and 177-00 white powder coated (RAL 9010)**

**Model 117-01 and 177-01 chromium-plated**

**Model 117-03 and 177-03 white aluminium powder coated (RAL 9006)**

**Model 117-04 and 177-04 grey aluminium powder coated (RAL 9007)**

German  
manufacture

All pictures:  
Model 177-04  
in basic equipment



### Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Possibility to pin on a paper roll holder on both traverse sides
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 2-sectioned:

- Head part upwards (+30°) adjustable by 2 metal rachets
- Length head part: 800 mm
- Cut-out for echocardiography and upholstery part to close, in patient supine position left, measures: L. 300 mm x W. 250 mm
- Fixed foot part
- Length foot part: 1150 mm
- Upholstery thickness: 80 mm

### Technical data:

- Length of lying surface: 1950 mm
- Width of lying surface: Model 117-...: 700 mm | Model 177-...: 800 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

### Mobility



Art. no. 1120-10  
Movable on 4 castors,  
Ø 100 mm, each with total  
lock




Art. no. 1140-10  
Movable on 4 castors,  
Ø 125 mm, each with  
total lock



Art. no. 1180N  
Movable on twin castors,  
Ø 100 mm with central lock  
(not in case of chromium-  
plated frame)



Art. no. 1185  
Movable on twin castors,  
Ø 125 mm with central lock  
(not in case of chromium-plated  
frame)

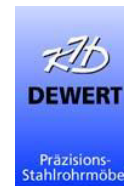

### Head part options



Art. no. 1655-10  
Head part stepless adjustable  
by gas spring, adjustment  
range: 0° - +60°

### Paper roll holder



Art. no. 1210 / 1220  
Paper roll holder, for width  
650 mm / 800 mm

### Upholstery options

Art. no. 1770 / 1780  
Soft covering for pleasant resting  
comfort, for width 650 mm / 800 mm,  
thickness 80 mm

Art. no. 1772 / 1782  
Upholstery medium-soft (cold foam  
cover), for width 650 mm / 800 mm,  
thickness 70 mm

### Side guards



Art. no. 1080-50  
1 side guard foldable,  
chromium-plated, completely made of steel  
tube, length: 1060 mm, length on cut-out  
side: 820 mm



Art. no. 1040  
1 side guard tiltable, chromium-plated, com-  
pletely made of steel tube, length: 700 mm



Art. no. 1005  
1 side guard to fit on side rail, with  
attachment clamp, chromium-plated,  
length: 800 mm,

### Add-on elements



Art. no. 1410  
1 piece of side rail,  
length: 350 mm

Art. no. 1420  
1 piece of side rail,  
length: 1300 mm



Art. no. 1430  
1 attachment clamp for  
side rail, suitable for pipe:  
Ø 18 mm



Art. no. 1440  
1 arm rest for side rail,  
without attachment clamp



Art. no. 1450  
Infusion pole, without  
attachment clamp



Art. no. 1850  
Storage basket with mounting,  
measures:  
L. 530 x W. 370 x H. 80 mm




# Tables for echocardiography

**2-sectioned table for echocardiography with electric motorized height adjustment ... for examination, treatment ...**

**Model 2550-00XLE and 2551-00XLE white powder coated (RAL 9010)**

**Model 2550-03XLE and 2551-03XLE white aluminium powder coated (RAL 9006)**

**Model 2550-04XLE and 2551-04XLE grey aluminium powder coated (RAL 9007)**

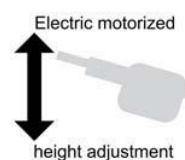

2 upper pictures:  
Model 2550-00XLE  
with options: Movable on twin  
castors with central  
lock (1180N),  
paper roll holder (1270K)



## Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- In head part: Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Plug connector for potential equalization in acc. to DIN 42801
- Dewert-hygiene-equipment



Model 2550-00XLE with options:  
Movable on comfort twin castors, each with lock  
(1120N-10), storage basket (1850), electric  
motorized head part adjustment (1665-10)

## Upholstery:

Lying surface 2-sectioned:

- Head part upwards (+25°) adjustable by 2 metal rachets
- Length head part: 800 mm
- In head part: Cut-out for echocardiography and upholstery part to close, in patient supine position left, measures: L. 300 mm x W. 250 mm
- Fixed foot part
- Length foot part: 1150 mm
- Upholstery thickness: 60 mm



## Technical data:


- Electronical height adjustment: 500 mm to 940 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: Model 2550-..XLE: 700 mm | Model 2551-..XLE: 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



## Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

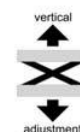



## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 1120N-10  
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm




Art. no. 1180N  
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185  
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N  
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering

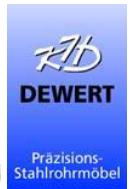

### Additional motor control options



Art. no. 1920E  
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E  
Foot switch (low voltage)



Art. no. 1940E  
1 foot switch fixation mounted at underframe

### Additional safety options



Art. no. 1950E  
Foot switch rail, lateral, on both sides, security adjustment



Art. no. B1997E  
Foot switch cover, made of brushed VA-steel



### Wall protection bumpers

Art. no. 1170  
1 wall protection bumper Ø 150 mm

### Head part options



Art. no. 1665-10  
Electric motorized head part adjustment to +60°



Art. no. 1655-10  
Head part stepless adjustable to +50° by gas spring

### Paper roll holder



Art. no. 1270K / 1280K  
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F  
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

### Upholstery options

Art. no. 1770 / 1780  
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782  
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

### Add-on elements



Art. no. 1410  
1 piece of side rail, length: 350 mm



Art. no. 1850  
Storage basket

### Side guards



Art. no. 1080-50  
1 side guard foldable, chromium-plated, completely made of steel tube, length: 1060 mm, length on cut-out side: 820 mm



Art. no. 1040  
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005  
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm,

# Tables for echocardiography

## 2-sectioned table for echocardiography with hydraulic height adjustment ... for examination, treatment ...


Model 2550-00XL/H and 2551-00XL/H white powder coated (RAL 9010)

Model 2550-03XL/H and 2551-03XL/H white aluminium powder coated (RAL 9006)

Model 2550-04XL/H and 2551-04XL/H grey aluminium powder coated (RAL 9007)



All pictures:  
Model 2550-00XL/H  
with options: Movable on  
twin castors with central  
lock (1180N), paper roll  
holder (1270K), head part  
adjustment by gas spring  
(1655-10)

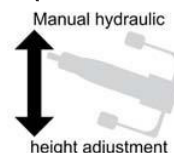

### Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- In head part: Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 2-sectioned:

- Head part upwards (+25°) adjustable by 2 metal rachets
- Length head part: 800 mm
- In head part: Cut-out for echocardiography and upholstery part to close, in patient supine position left, measures: L. 300 mm x W. 250 mm
- Fixed foot part
- Length foot part: 1150 mm
- Upholstery thickness: 60 mm



### Technical data:


- Manual foot operated height adjustment: 520 mm to 940 mm by hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: Model 2550-...XL/H: 700 mm | Model 2551-...XL/H: 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

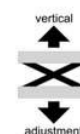

## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 1120N-10  
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N  
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185  
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N  
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



### Head part options



Art. no. 1655-10  
Head part stepless adjustable to +50° by gas spring



Art. no. 1270K / 1280K  
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F  
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

### Add-on elements



Art. no. 1410  
1 piece of side rail, length: 350 mm



Art. no. 1850  
Storage basket

### Side guards



Art. no. 1080-50  
1 side guard foldable, chromium-plated, completely made of steel tube, length: 1060 mm, length on cut-out side: 820 mm



Art. no. 1040  
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005  
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm,

### Wall protection bumpers



Art. no. 1170  
1 wall protection bumper Ø 150 mm

### Upholstery options

Art. no. 1770 / 1780  
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782  
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm




# Tables for echocardiography

**3-sectioned table for echocardiography with electric motorized height adjustment ... for examination, treatment ...**

**Model 2552-00XLE white powder coated (RAL 9010)**

**Model 2552-03XLE white aluminium powder coated (RAL 9006)**

**Model 2552-04XLE grey aluminium powder coated (RAL 9007)**

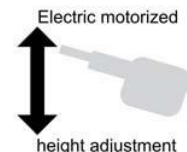

**NEW**

German  
manufacture

without  
picture



Model 2552-00XL/H with options:  
Movable on twin castors with central lock (1180N),  
head part adjustment by gas spring (1655-10),  
paper roll holder (1270K)



Frame in case of electric motorized height adjustment

## Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- In head part: Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Plug connector for potential equalization in acc. to DIN 42801
- Dewert-hygiene-equipment



## Upholstery:

Lying surface 3-sectioned:

- Head part upwards (+25°) adjustable by 2 metal rachets
- Length head part: 800 mm
- In head part: Cut-out for echocardiography and upholstery part to close, in patient supine position left, measures: L. 300 mm x W. 250 mm
- Fixed middle part, length: 300 mm
- Foot part stepless upwards (+30°) adjustable by gas spring
- Length foot part: 850 mm
- Upholstery thickness: 60 mm

## Technical data:


- Electronical height adjustment: 540 mm to 980 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm, **optional** 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



## Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

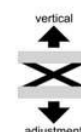

## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 1120N-10  
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N  
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185  
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N  
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



### Additional motor control options



Art. no. 1920E  
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E  
Foot switch (low voltage)



Art. no. 1940E  
1 foot switch fixation mounted at underframe

### Additional safety options



Art. no. 1950E  
Foot switch rail, lateral, on both sides, security adjustment



Art. no. B1997E  
Foot switch cover, made of brushed VA-steel



### Wall protection bumpers

Art. no. 1170  
1 wall protection bumper Ø 150 mm

### Head part options



Art. no. 1665-10 electric motorized head part adjustment to +60°



Art. no. 1655-10  
Head part stepless adjustable to +50° by gas spring

### Paper roll holder



Art. no. 1270K / 1280K  
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F  
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

### Upholstery options

Art. no. 2552-800  
Width of upholstery: 800 mm

Art. no. 1770 / 1780  
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782  
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

### Add-on elements



Art. no. 1410  
1 piece of side rail, length: 350 mm



Art. no. 1850  
Storage basket

### Side guards



Art. no. 1080-50  
1 side guard foldable, chromium-plated, completely made of steel tube, length: 1060 mm, length on cut-out side: 820 mm



Art. no. 1040  
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005  
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm,


# Tables for echocardiography

## 3-sectioned table for echocardiography with hydraulic height adjustment ... for examination, treatment ...

Model 2552-00XL/H white powder coated (RAL 9010)

Model 2552-03XL/H white aluminium powder coated (RAL 9006)

Model 2552-04XL/H grey aluminium powder coated (RAL 9007)

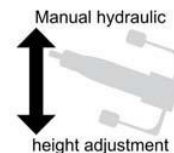

**NEU**

German  
manufacture

All pictures:  
Model 2552-00XL/H  
with options: Movable on  
twin castors with central  
lock (1180N), head part  
adjustment by gas spring  
(1655-10), paper roll holder  
(1270K)



### Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- In head part: Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 3-sectioned:

- Head part upwards (+25°) adjustable by 2 metal rachets
- Length head part: 800 mm
- In head part: Cut-out for echocardiography and upholstery part to close, in patient supine position left, measures: L. 300 mm x W. 250 mm
- Fixed middle part, length: 300 mm
- Foot part stepless upwards (+30°) adjustable by gas spring
- Length foot part: 850 mm
- Upholstery thickness: 60 mm



### Technical data:


- Manual foot operated height adjustment: 560 mm to 980 mm by hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm, **optional** 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

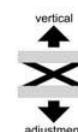



## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 1120N-10  
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N  
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185  
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N  
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



### Head part options



Art. no. 1655-10  
Head part stepless adjustable to +50° by gas spring



Art. no. 1270K / 1280K  
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F  
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

### Add-on elements



Art. no. 1410  
1 piece of side rail, length: 350 mm



Art. no. 1850  
Storage basket

### Side guards



Art. no. 1080-50  
1 side guard foldable, chromium-plated, completely made of steel tube, length: 1060 mm, length on cut-out side: 820 mm



Art. no. 1040  
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005  
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm,

### Wall protection bumpers



Art. no. 1170  
1 wall protection bumper Ø 150 mm

### Upholstery options

Art. no. 2552-800  
Width of upholstery: 800 mm

Art. no. 1770 / 1780  
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782  
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

# Tables for echocardiography


## 3-sectioned table for echocardiography with electric motorized height adjustment

... for examination, treatment ...

**Model 2690-00XLE** white powder coated (RAL 9010)

**Model 2690-03XLE** white aluminium powder coated (RAL 9006)

**Model 2690-04XLE** grey aluminium powder coated (RAL 9007)

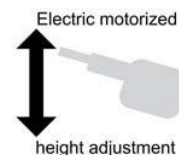

All pictures:  
Model 2690-04XLE  
with option:  
Movable on twin castors with  
central lock 1187); with  
special feature: foot bar



### Basic equipment:

- Table movable on twin castors, Ø 75 mm, each with lock
- In middle part: Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Middle part stepless positive and negative adjustable by electromotive drive, control by hand switch
- Possible positions including: Roof position, Trendelenburg, Reverse Trendelenburg
- Plug connector for potential equalization in acc. to DIN 42801
- Electromotive drive control of height adjustment by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 3-sectioned:

- Head part stepless positive (+45°) and negative (-30°) adjustable by gas spring
- Length head part: 400 mm, with nose opening and upholstery part to close
- Middle part (length: 580 mm) stepless positive (+15°) and negative (-15°) adjustable by electromotive drive
- Middle part with cut-out for echocardiography and upholstery part to close, measures: L. 300 mm x W. 250 mm
- Foot part stepless positive (+62°) and negative (-12°) adjustable by gas spring
- Length foot part: 970 mm
- Upholstery thickness: 60 mm



### Technical data:


- Electronical height adjustment: 630 mm to 1070 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 220 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

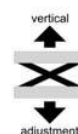

Available as an option (see pricelist for additional options)

### Mobility



Art. no. 1182  
Without basic mobility



Art. no. 1187  
Movable on twin castors with central lock,  
Ø 100 mm, lower underframe,  
table height 540 mm – 980 mm

### Additional motor control options



Art. no. 1920E  
Hand switch (low voltage),  
with magnet on the back to fasten



Art. no. 1910E  
Foot switch (low voltage)



Art. no. 1940E  
1 foot switch fixation mounted at  
underframe

### Additional safety options



Art. no. B1997E  
Foot switch cover,  
made of brushed VA-steel

### Nose opening



Art. no. 1615-EF  
Without nose opening

### Wall protection bumpers



Art. no. 1175  
4 wall protection bumpers Ø 150 mm,  
including mounting, for width 700 mm

### Paper roll holder



Art. no. 1270F/K  
Paper roll holder, mounted under head  
part (not protruding),  
for width 700 mm



Art. no. 1270F  
Paper roll holder, mounted at foot part  
(protruding: 70 mm),  
for width 700 mm

### Add-on elements



Art. no. 1410  
1 piece of side rail, length: 350 mm



Art. no. 1005  
1 side guard to fit on side rail, with  
attachment clamp, chromium-plated,  
length: 800 mm,

### Upholstery options

Art. no. 1772  
Upholstery medium-soft (cold foam  
cover), for width 700 mm,  
thickness 70 mm



# Tables for echocardiography


## 3-sectioned table for echocardiography with electric motorized height adjustment

... for examination, treatment ...

**Model 2695-00XLE** white powder coated (RAL 9010)

**Model 2695-03XLE** white aluminium powder coated (RAL 9006)

**Model 2695-04XLE** grey aluminium powder coated (RAL 9007)

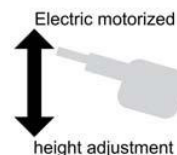

German  
manufacture

All pictures:  
Model 2695-00XLE  
with option: Paper roll  
holder (1270F/K)



### Basic equipment:

- Table movable on twin castors, Ø 75 mm, each with lock
- In middle part: 2 cut-outs for echocardiography and upholstery parts to close
- Middle part stepless positive and negative adjustable by electromotive drive, control by hand switch
- Possible positions including: Roof position, Trendelenburg, Reverse Trendelenburg
- Plug connector for potential equalization in acc. to DIN 42801
- Electromotive drive control of height adjustment by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 3-sectioned:

- Head part stepless positive (+45°) and negative (-30°) adjustable by gas spring
- Length head part: 400 mm, with nose opening and upholstery part to close
- Middle part (length: 580 mm) stepless positive (+20°) and negative (-10°) adjustable by electromotive drive
- Middle part with cut-out for echocardiography and upholstery part to close, measures: L. 300 mm x W. 250 mm
- Foot part stepless positive (+60°) and negative (-12°) adjustable by gas spring
- Length foot part: 970 mm
- Upholstery thickness: 60 mm



### Technical data:

- Electronical height adjustment: 630 mm to 1070 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 220 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

### Mobility



Art. no. 1182  
Without basic mobility



Art. no. 1187  
Movable on twin castors with central lock,  
Ø 100 mm, lower underframe,  
table height 540 mm – 980 mm



### Additional motor control options



Art. no. 1920E  
Hand switch (low voltage),  
with magnet on the back to fasten



Art. no. 1910E  
Foot switch (low voltage)



Art. no. 1940E  
1 foot switch fixation mounted at  
underframe

### Additional safety options



Art. no. B1997E  
Foot switch cover,  
made of brushed VA-steel



Art. no. 1615-EF  
Without nose opening

### Paper roll holder



Art. no. 1270F/K  
Paper roll holder, mounted under head  
part (not protruding),  
for width 700 mm



Art. no. 1270F  
Paper roll holder, mounted at foot part  
(protruding: 70 mm),  
for width 700 mm

### Wall protection bumpers



Art. no. 1175  
4 wall protection bumpers Ø 150 mm,  
including mounting, for width 700 mm

### Add-on elements



Art. no. 1410  
1 piece of side rail, length: 350 mm



Art. no. 1005  
1 side guard to fit on side rail, with  
attachment clamp, chromium-plated,  
length: 800 mm,

### Upholstery options

Art. no. 1772  
Upholstery medium-soft (cold foam cover),  
for width 700 mm, thickness 70 mm




# Tilt table for echokardiography

**Tilt table with electric motorized height adjustment and tilt function**  
... for examination, treatment ...

**Model 2095-00E white powder coated (RAL 9010)**

**Model 2095-03E white aluminium powder coated (RAL 9006)**

**Model 2095-04E grey aluminium powder coated (RAL 9007)**

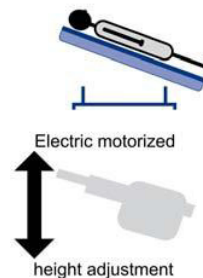

Pictures:

**Model 2090-04E with options**

(1 pair of big grip handles (2090-038), head part manual upwards adjustable (2090-058), 1 pair of belt guides (2090-036), 1 safety belt with automatic lock (2090-0351), paper roll holder (1272))

without picture

German manufacture



Model 2095-04E differs from model 2090-04E (see picture):

- Cut-out for echocardiography, conical shape prevents the unintended slipping out of upholstery part
- Head part adjustment by gas spring

## Basic equipment:

- Lying surface 2-sectioned, cut-out for echocardiography with upholstery part to close in head part
- Operation of electric motorized height adjustment by foot switch, operation of electric motorized tilt function by hand switch (independent from each other)
- Integrated safety function for automatic blocking of electromotive drive
- Plug connector for potential equalization in acc. to DIN 42801
- Movable on 4 twin castors Ø 75 mm, each with lock
- Foot plate to stand on (upholstered), removable
- Side rails along lying surface (left and right side) for special tilt table equipment
- Dewert-hygiene-equipment



## Upholstery:

Lying surface 2-sectioned:

- Head part stepless adjustable from 0° to +45° by gas spring, length head part: 800 mm
- In head part: Cut-out for echocardiography, measures: Length 300 mm x Width 250 mm
- Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Fixed foot part, length foot part: 1130 mm
- Thickness of upholstery: 80 mm

## Technical data:

- Electric motorized height adjustment: 620 mm to 930 mm in 34 sec. by AC-Motor, pressure capacity of 10.000 N
- Electric motorized tilt function from 0° to +85° in 22 sec. by AC-Motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm | Width overall: 850 mm
- Length of lying surface: 1930 mm | Length overall: 1960 mm
- Load capacity: 180 kg; according to IEC 60601 tested with fourfold static load



## Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology





Available as an option (see pricelist for additional options)

### Additional motor control options



Art. no. 1920E  
Hand switch (low voltage),  
with magnet on the back to  
fasten



Art. no. 1910E  
1 foot switch (low voltage)

### Grip handles



Art. no. 2090-028  
1 pair grip handles, small



Art. no. 2090-038  
1 pair grip handles, big

### Belt elements



Art. no. 2090-034  
Upholstery for strap



Art. no. 2090-035  
1 safety belt with velcro fastener



Art. no. 2090-0351  
1 safety belt with automatic  
lock



Art. no. 2090-036  
1 pair belt guides

### Arm rest



Art. no. 2090-0143  
Attachment clamp to fix  
arm rest at side rail

Art. no. 2090-059  
Arm rest

### Paper roll holder



Art. no. 1272  
Paper roll holder, mounted under head  
part ( not protruding)

### Special equipment



Art. no. 2090-090  
Angle gauge, on both sides  
of the table

K.H. Dewert GmbH  
Vollmestraße 7  
33649 Bielefeld  
Germany

Post office 14 40 28  
33634 Bielefeld  
Germany

Phone + 49 (0) 5 21 / 4 00 27-0  
Fax + 49 (0) 5 21 / 4 00 27-27

Internet: [www.khdewert.de](http://www.khdewert.de)  
E-Mail: [info@khdewert.de](mailto:info@khdewert.de)  
YouTube: K.H. Dewert

## TABLES FOR ECHOCARDIOGRAPHY FOR EXAMINATION AND TREATMENT

TABLES FOR SURGERY AND HOSPITAL

