


THERAPY TABLES

2021.1



# SERIES 2800XL AND 4800E

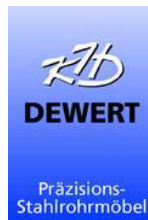

INCLINATION  
ADJUSTMENT




ROOF POSITION



GYMNASTIC TABLES  
(BOBATH & VOJTA)



THERAPY TABLES

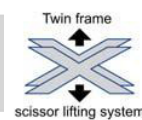



# Therapy tables

## Content

Good reasons for a therapy table of series 2800XL and 4800E	3
---	---

### Therapy tables - Series 2800XL

#### 3-sectioned therapy tables:

• Model 2815XLE	6
• Model 2815XL/H	8
• Model 2820XLE	10
• Model 2820XL/H	12



#### 3-sectioned therapy tables, all upholstery sections adjustable:

• Model 2830XLE	14
• Model 2830XL/H	16
• Model 2890XLE	18
• Model 2890XL/H	20



#### 4-sectioned therapy tables, all upholstery sections adjustable:

• Model 2895XLE	22
• Model 2895XL/H	24
• Model 2850XLE	26
• Model 2850XL/H	28
• Model 2860XLE	30
• Model 2860XL/H	32



#### 6-sectioned therapy tables, all upholstery sections adjustable:

• Model 2880XLE	34
• Model 2880XL/H	36



### Therapy tables - Series 4800E

• Model 4815E	38
• Model 4840E	40

#### Gymnastic tables for therapy acc. to Bobath and Vojta

• Series 2800XLE	42
• Series 4800E	44




Additional equipment	46
----------------------	----



All therapy tables of series 2800XL and 4800E with the exceptions of gymnastic tables have a nose opening with upholstery part to close in head part as a standard.


# Therapy tables

## Good reasons for a therapy table of series 2800XL and 4800E

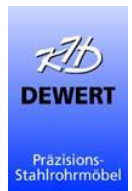


### Electric motorized height adjustment Series 2800XL und 4800E

- Electronical height adjustment by AC-motor with a pressure capacity of 10000 N
- Constant adjustment speed even in case of high load
- EPR-safety power cable with strain relief, length: 4 m
- **Integrated safety devices (according to the newest advice of BfArM)**
  - Automatic blocking of electromotive drive
  - Active freewheeling downwards, automatic declutch in case of resistance
- High load capacity
- Series 4800E additionally: High adjustment range

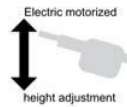

### Hydraulic height adjustment Series 2800XL

- Manual foot operated height adjustment by a hydraulic pump with a pressure capacity of 10000 N
- Pedals for foot-controlled height positioning on both sides of the table
- Load capacity: 250 kg

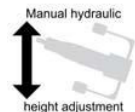

### Dewert-hygiene-equipment Enables optimal cleaning and disinfection




- Hinge lid made of identical upholstery covering
- Artificial leather underneath the upholstery
- Vent holes underneath the upholstery, closed by special caps with valve function
- Easy hygiene because of an open construction: good accessibility of components, no plastic coverings which are breakable in the ageing process
- Option: Innovative artificial leather Skai® Toronto EN with staynu






# Therapy tables


## Good reasons for a therapy table of series 2800XL and 4800E




**Durable** ✓




**Robust** ✓




**Easy-care** ✓

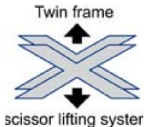


**Twin frame scissor lifting system** ✓




**High load capacity** ✓



**Proven Dewert-quality**

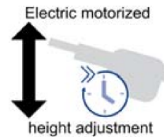

**High adjustment speed** ✓



**Universal use** ✓

### Technology:




- Extremely stable twin frame scissor lifting system, with high torsional stiffness due to thick-walled square tubes with big dimensions, also for under and upper frame and scissor



- High stability due to massive 20 mm axes
- Widely-spaced feet and a high deadweight guarantee stability
- Robust frame construction which is proven in demanding everyday life of hospitals for years



- Abrasion-resistant guide system at under and upper frame due to the use of roller bearings

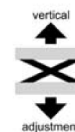

### Upholstery:

- Thickness of upholstery: 60 mm
- Multiple glued, high-strength laminated wood panel (thickness: 18 mm) no breaking chipboards in case of pressure on the edges
- Two-layer homogenous polyurethane foam, density 46 and 75, Ökotex-certified
- 2 different variants for artificial leather eligible for upholstery covering:
  - Skai® Pandoria Plus with velvet-smooth surface
  - Skai® Toronto EN with innovative staynu-technology



Cleaning kit for Skai products  
(Art. no. 050-1000)

Consists of:

- Spray bottle with 200 ml special Skai cleansing agent
- 1 microfibre care cloth 40 x 40 cm

# Therapy tables

## Good reasons for a therapy table of series 2800XL and 4800E




### ... with 3-sectioned lying surface:

- Adjustable head part either positive and negative or only positive adjustable with tubular steel frame for high torsion rigidity
- Pressure stable gas springs ensure working at head part and foot part **without vibrations**
- Operating lever to manipulate the gas spring easily accessible
- **Variable length** due to foldable head part (partly: can be hinged vertical)



### ... with 3-, 4- or 6-sectioned lying surface, all upholstery sections adjustable:

- Head part positive and negative adjustable
- Adjustment of upholstery sections either by pressure stable gas springs or partly by a self-fixing system of 2 guide rods, without clearance (depending on model)
- Operating levers to manipulate the gas spring easily accessible
- Middle part adjustable to **roof position**, **Trendelenburg position** possible
- **Models 2890 and 2895:**
  - Multifunctional and independent adjustment of upholstery sections:  
**Roof position, Trendelenburg, reverse Trendelenburg**

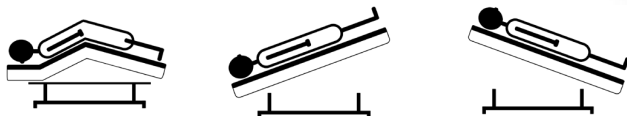

### ... with 1-sectioned lying surface:

- Gymnastic tables for therapy according to Bobath and Vojta
- Width of lying surface eligible:  
100 cm, 110 cm and 120 cm




# Series 2800XL

## 3-sectioned therapy table with electric motorized height adjustment

Model 2815-00XLE white powder coated (RAL 9010)

Model 2815-03XLE white aluminium powder coated (RAL 9006)

Model 2815-04XLE grey aluminium powder coated (RAL 9007)

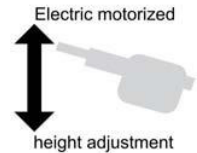

without  
picture



Model 2815-00XL/H (hydraulic height adjustment) with option:  
Movable on twin castors with central lock (046)



### Basic equipment:

- Lying surface 3-sectioned
- Foot part stepless upwards (+75°) adjustable to sitting position
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Under part in case  
of electric motorized  
height adjustment



### Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+45°) and downwards (-35°) adjustable by gas spring
- Length head part: 540 mm, with nose opening and upholstery part to close
- Fixed middle part
- Length middle part: 530 mm
- Foot part stepless upwards (+75°) adjustable to sitting position by gas spring
- Length foot part: 880 mm
- Thickness of upholstery: 60 mm



### Technical data:


- Electronical height adjustment: 510 mm to 950 mm in 22 sec.  
by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

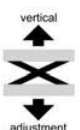



## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 040N  
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. B1185  
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



### Additional motor control options



Art. no. 090E  
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 091E  
Foot switch (low voltage)



Art. no. 092E  
1 foot switch fixation mounted at underframe



Art. no. 094E  
Foot switch rail, lateral, on both sides, security adjustment



Art. no. 095E  
Foot switch rail, surround, security adjustment

### Additional safety options



Art. no. B1997E  
Foot switch cover, made of brushed VA-steel

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery part to close

### Head part options



Art. no. 056  
Head part 3-sectioned, arm rests downwards adjustable from -25° to +45° and removable (art. no. 010K not possible)  
measures:  
Head part: L. 540 mm x W. 370 mm  
Arm rests: L. 500 mm x W. 160 mm

### Upholstery options

Art. no. 007  
Width of upholstery 800 mm

Art. no. B1770  
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. B1772  
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

### Paper roll holder



Art. no. 010K  
Paper roll holder, mounted under head part (not protruding), for width 700 mm, not in case of 3-sectioned head part (Art. no. 056)



Art. no. 010K/F  
Paper roll holder, mounted at foot part (not protruding), for width 700 mm

### Add-on elements



Art. no. B1410  
1 piece of side rail, length: 350 mm



Art. no. B1430  
1 attachment clamp for side rail, suitable for pipe: Ø 18 mm

### Arm rests

without picture

Art. no. 088  
1 arm rest, fixed, swivelling and removeable, additional revolving and tiltable



Art. no. B1440  
1 arm rest for side rail, without attachment clamp, swivelling and removeable, additional revolving and tiltable


# Series 2800XL

## 3-sectioned therapy table with hydraulic height adjustment

**Model 2815-00XL/H white powder coated (RAL 9010)**

**Model 2815-03XL/H white aluminium powder coated (RAL 9006)**

**Model 2815-04XL/H grey aluminium powder coated (RAL 9007)**

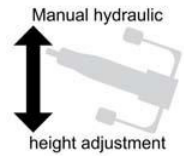

All pictures:  
Model 2815-00XL/H  
with options:  
Movable on twin castors with  
central lock (046)



### Basic equipment:

- Lying surface 3-sectioned
- Foot part stepless upwards (+75°) adjustable to sitting position
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+45°) and downwards (-35°) adjustable by gas spring
- Length head part: 540 mm, with nose opening and upholstery part to close
- Fixed middle part
- Length middle part: 530 mm
- Foot part stepless upwards (+75°) adjustable to sitting position by gas spring
- Length foot part: 880 mm
- Thickness of upholstery: 60 mm



### Technical data:


- Manual foot operated height adjustment: 530 mm to 950 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

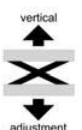



## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 040N  
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. B1185  
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



### Paper roll holder



Art. no. 010K  
Paper roll holder, mounted under head part (not protruding), for width 700 mm, not in case of 3-sectioned head part (Art. no. 056)



Art. no. 010K/F  
Paper roll holder, mounted at foot part (not protruding), for width 700 mm

### Upholstery options

Art. no. 007  
Width of upholstery 800 mm

Art. no. B1770  
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. B1772  
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

### Head part options



Art. no. 056  
Head part 3-sectioned, arm rests downwards adjustable from -25° to +45° and removable (art. no. 010K not possible)  
measures:  
Head part: L. 540 mm x W. 370 mm  
Arm rests: L. 500 mm x W. 160 mm

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery part to close

### Add-on elements




Art. no. B1410  
1 piece of side rail, length: 350 mm



Art. no. B1430  
1 attachment clamp for side rail, suitable for pipe: Ø 18 mm

### Arm rests

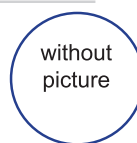

Art. no. 088  
1 arm rest, fixed, swivelling and removeable, additional revolving and tiltable



Art. no. B1440  
1 arm rest for side rail, without attachment clamp, swivelling and removeable, additional revolving and tiltable


# Series 2800XL

## 3-sectioned therapy table with electric motorized height adjustment

Model 2820-00XLE white powder coated (RAL 9010)

Model 2820-03XLE white aluminium powder coated (RAL 9006)

Model 2820-04XLE grey aluminium powder coated (RAL 9007)

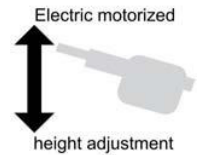

All pictures:  
Model 2820-04XLE  
with options:  
Movable on twin castors with  
central lock (046), side rail  
(B1410)



### Basic equipment:

- Lying surface 3-sectioned
- Head part can be hinged vertical
- Foot part stepless upwards adjustable to sitting position (+80°)
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length head part: 450 mm, with nose opening and upholstery part to close
- Fixed middle part
- Length middle part: 530 mm
- Foot part stepless upwards (+80°) adjustable to sitting position by gas spring
- Length foot part: 1020 mm
- Thickness of upholstery: 60 mm



### Technical data:

- Electronical height adjustment: 510 mm to 950 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 2000 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 040N  
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. B1185  
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



### Additional motor control options



Art. no. 090E  
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 091E  
Foot switch (low voltage)



Art. no. 092E  
1 foot switch fixation mounted at underframe

### Additional safety options



Art. no. 094E  
Foot switch rail, lateral, on both sides, security adjustment



Art. no. B1997E  
Foot switch cover, made of brushed VA-steel

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery part to close

### Paper roll holder



Art. no. 010K  
Paper roll holder, mounted under head part (not protruding), for width 700 mm, not in case of 3-sectioned head part (Art. no. 056)



Art. no. 010K/F  
Paper roll holder, mounted at foot part (not protruding), for width 700 mm

### Upholstery options

Art. no. 007  
Width of upholstery 800 mm

Art. no. B1770  
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. B1772  
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

### Add-on elements



Art. no. B1410  
1 piece of side rail, length: 350 mm



Art. no. B1430  
1 attachment clamp for side rail, suitable for pipe: Ø 18 mm

### Arm rests



Art. no. B1440  
1 arm rest for side rail, without attachment clamp, swivelling and removeable, additional revolving and tiltable




# Series 2800XL

## 3-sectioned therapy table with hydraulic height adjustment

**Model 2820-00XL/H white powder coated (RAL 9010)**

**Model 2820-03XL/H white aluminium powder coated (RAL 9006)**

**Model 2820-04XL/H grey aluminium powder coated (RAL 9007)**

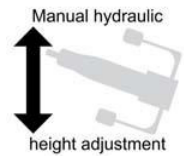

All pictures:  
Model 2820-04XL/H  
with options:  
Movable on twin castors with  
central lock (046), paper roll  
holder (010K/F)



### Basic equipment:

- Lying surface 3-sectioned
- Head part can be hinged vertical
- Foot part stepless upwards adjustable to sitting position (+80°)
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length head part: 450 mm, with nose opening and upholstery part to close
- Fixed middle part
- Length middle part: 530 mm
- Foot part stepless upwards (+80°) adjustable to sitting position by gas spring
- Length foot part: 1020 mm
- Thickness of upholstery: 60 mm



### Technical data:


- Manual foot operated height adjustment: 530 mm to 950 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 2000 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

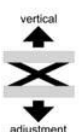

Available as an option (see pricelist for additional options)

### Mobility



Art. no. 040N  
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. B1185  
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



### Paper roll holder



Art. no. 010K  
Paper roll holder, mounted under head part (not protruding), for width 700 mm



Art. no. 010K/F  
Paper roll holder, mounted at foot part (not protruding), for width 700 mm

### Upholstery options

Art. no. 007  
Width of upholstery 800 mm

Art. no. B1770  
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. B1772  
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

### Add-on elements



Art. no. B1410  
1 piece of side rail, length: 350 mm



Art. no. B1430  
1 attachment clamp for side rail, suitable for pipe: Ø 18 mm

### Arm rests



Art. no. B1440  
1 arm rest for side rail, without attachment clamp, swivelling and removeable, additional revolving and tiltable

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery part to close


# Series 2800XL

## 3-sectioned therapy table with electric motorized height adjustment

Model 2830-00XLE white powder coated (RAL 9010)

Model 2830-03XLE white aluminium powder coated (RAL 9006)

Model 2830-04XLE grey aluminium powder coated (RAL 9007)

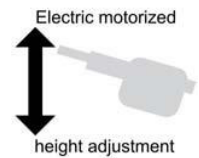

All pictures: Model 2830-00XLE with options:  
Wheel-lifting/fixing mechanism (045N), surrounding foot switch rail (095E), 3-sectioned head part (056)



### Basic equipment:

- Lying surface 3-sectioned
- Middle part stepless adjustable to roof position, Trendelenburg position possible
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part stepless adjustable (+70°) by a self-fixing system of 2 guide rods, without clearance
- Length foot part: 940 mm
- Thickness of upholstery: 60 mm



### Technical data:


- Electronical height adjustment: 510 mm to 950 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

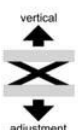



## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 040N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock,  
basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm with  
central lock, basic height +20 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering



### Additional motor control options



Art. no. 090E  
Hand switch (low voltage),  
with magnet on the back to fasten



Art. no. 091E  
Foot switch (low voltage)



Art. no. 092E  
1 foot switch fixation mounted at  
underframe



Art. no. 094E  
Foot switch rail, lateral, on both sides,  
security adjustment



Art. no. 095E  
Foot switch rail, surround,  
security adjustment



Art. no. 065E  
Electric motorized adjustment to roof  
position

### Additional safety options



Art. no. B1997E  
Foot switch cover,  
made of brushed VA-steel

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery  
part to close

### Head part options



Art. no. 056  
Head part 3-sectioned, arm rests down-  
wards adjustable from -25° to +45° and  
removable (art. no. 010K not possible)  
measures:  
Head part: L. 540 mm x W. 370 mm  
Arm rests: L. 500 mm x W. 160 mm

### Paper roll holder



Art. no. 010K  
Paper roll holder, mounted under head  
part (not protruding), for width 700 mm,  
not in case of 3-sectioned head part  
(Art. no. 056)



Art. no. 010K/F  
Paper roll holder, mounted at foot part  
(not protruding), for width 700 mm

### Arm rests

without  
picture

Art. no. 088  
1 arm rest, fixed, swivelling and  
removeable, additional revolving and  
tiltable

### Upholstery options

Art. no. 007  
Width of upholstery 800 mm

Art. no. B1772  
Upholstery medium-soft (cold foam cover),  
for width 700 mm, thickness 70 mm


# Series 2800XL

## 3-sectioned therapy table with hydraulic height adjustment

Model 2830-00XL/H white powder coated (RAL 9010)

Model 2830-03XL/H white aluminium powder coated (RAL 9006)

Model 2830-04XL/H grey aluminium powder coated (RAL 9007)

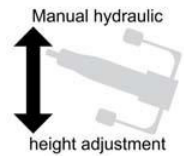

All pictures:  
Model 2830-00XL/H  
with options:  
Wheel-lifting/fixing mechanism  
(045N), paper roll holder (010K)



### Basic equipment:

- Lying surface 3-sectioned
- Middle part stepless adjustable to roof position, Trendelenburg position possible
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part stepless adjustable (+70°) by a self-fixing system of 2 guide rods, without clearance
- Length foot part: 940 mm
- Thickness of upholstery: 60 mm



### Technical data:


- Manual foot operated height adjustment: 530 mm to 950 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

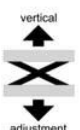

Available as an option (see pricelist for additional options)

### Mobility



Art. no. 040N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock,  
basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm with  
central lock, basic height +20 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering



### Paper roll holder



Art. no. 010K  
Paper roll holder, mounted under head  
part (not protruding), for width 700 mm,  
not in case of 3-sectioned head part  
(Art. no. 056)



Art. no. 010K/F  
Paper roll holder, mounted at foot part  
(not protruding), for width 700 mm




Art. no. 056  
Head part 3-sectioned, arm rests  
downwards adjustable from -25° to  
+45° and removable  
(art. no. 010K not possible)  
measures:  
Head part: L. 540 mm x W. 370 mm  
Arm rests: L. 500 mm x W. 160 mm

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery  
part to close

### Arm rests

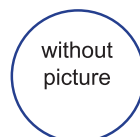

Art. no. 088  
1 arm rest, fixed, swivelling and  
removeable, additional revolving and  
tiltable

### Upholstery options

Art. no. 007  
Width of upholstery 800 mm

Art. no. B1772  
Upholstery medium-soft (cold foam  
cover), for width 700 mm,  
thickness 70 mm




# Series 2800XL

## 3-sectioned therapy table with electric motorized height adjustment

**Model 2890-00XLE white powder coated (RAL 9010)**

**Model 2890-03XLE white aluminium powder coated (RAL 9006)**

**Model 2890-04XLE grey aluminium powder coated (RAL 9007)**

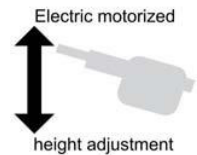

without picture



Under frame in case of electric motorized height adjustment, with options:  
Movable on twin castors with central lock (046)

Model 2890-00XL/H (hydraulic height adjustment)  
with options:  
Movable on twin castors with central lock (046),  
paper roll holder (010K/F)

### Basic equipment:

- Lying surface 3-sectioned
- Multifunctional and independent adjustment of upholstery sections
- Possibilities of positioning the surface parts: Roof shape, Trendelenburg, reverse Trendelenburg
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 420 mm; with nose opening and upholstery part to close
- Middle part stepless up- (+45°) and downwards (-15°) adjustable by gas spring
- Length middle part: 550 mm
- Foot part stepless up- (+50°) and downwards (-15°) adjustable by gas spring
- Length foot part: 980 mm
- Thickness of upholstery: 60 mm

### Technical data:


- Electronical height adjustment: 570 mm to 1010 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

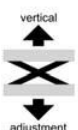

## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 077N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock



Art. no. 049  
Movable on twin castors, Ø 100 mm  
with central lock



Art. no. 044N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering



### Additional motor control options



Art. no. 090E  
Hand switch (low voltage),  
with magnet on the back to fasten



Art. no. 091E  
Foot switch (low voltage)



Art. no. 092E  
1 foot switch fixation mounted at  
underframe

### Additional safety options



Art. no. B1997E  
Foot switch cover,  
made of brushed VA-steel

### Paper roll holder



Art. no. 010F/K  
Paper roll holder, mounted under head  
part (not protruding), for width 700 mm



Art. no. 010F  
Paper roll holder, mounted at foot part  
(protruding: 70 mm), for width 700 mm

### Add-on elements



Art. no. B1410  
1 piece of side rail, length: 350 mm



Art. no. B1430  
1 attachment clamp for side rail,  
suitable for pipe: Ø 18 mm

### Arm rests



Art. no. B1440  
1 arm rest for side rail, without  
attachment clamp, swivelling and  
removeable, additional revolving and  
tiltable

### Side guards



Art. no. 1005  
1 side guard to fit on side rail, with  
attachment clamp, chromium-plated,  
length: 800 mm

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery  
part to close

### Upholstery options

Art. no. B1770  
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. B1772  
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm


# Series 2800XL

## 3-sectioned therapy table with hydraulic height adjustment

**Model 2890-00XL/H white powder coated (RAL 9010)**

**Model 2890-03XL/H white aluminium powder coated (RAL 9006)**

**Model 2890-04XL/H grey aluminium powder coated (RAL 9007)**

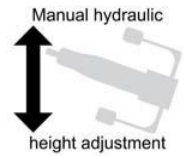

All pictures:  
Model 2890-00XL/H  
with options:  
Movable on twin castors with  
central lock (046), paper roll  
holder (010K/F)



### Basic equipment:

- Lying surface 3-sectioned
- Multifunctional and independent adjustment of upholstery sections
- Possibilities of positioning the surface parts: Roof shape, Trendelenburg, reverse Trendelenburg
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 420 mm; with nose opening and upholstery part to close
- Middle part stepless up- (+45°) and downwards (-15°) adjustable by gas spring
- Length middle part: 550 mm
- Foot part stepless up- (+50°) and downwards (-15°) adjustable by gas spring
- Length foot part: 980 mm
- Thickness of upholstery: 60 mm



### Technical data:


- Manual foot operated height adjustment: 600 mm to 1020 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

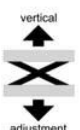



## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 077N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock



Art. no. 049  
Movable on twin castors, Ø 100 mm  
with central lock



Art. no. 044N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering



### Paper roll holder



Art. no. 010F/K  
Paper roll holder, mounted under head  
part (not protruding), for width 700 mm



Art. no. 010F  
Paper roll holder, mounted at foot part  
(protruding: 70 mm), for width 700 mm



Art. no. 1615-EF  
Loss of nose opening and upholstery  
part to close

### Add-on elements



Art. no. B1410  
1 piece of side rail, length: 350 mm



Art. no. B1430  
1 attachment clamp for side rail,  
suitable for pipe: Ø 18 mm

### Arm rests



Art. no. B1440  
1 arm rest for side rail, without  
attachment clamp, swivelling and  
removeable, additional revolving and  
tiltable

### Side guards



Art. no. 1005  
1 side guard to fit on side rail, with  
attachment clamp, chromium-plated,  
length: 800 mm

### Upholstery options

Art. no. B1770  
Soft covering for pleasant resting comfort,  
for width 700 mm, thickness 80 mm

Art. no. B1772  
Upholstery medium-soft (cold foam cover),  
for width 700 mm, thickness 70 mm


# Series 2800XL

## 4-sectioned therapy table with electric motorized height adjustment

**Model 2895-00XLE** white powder coated (RAL 9010)

**Model 2895-03XLE** white aluminium powder coated (RAL 9006)

**Model 2895-04XLE** grey aluminium powder coated (RAL 9007)

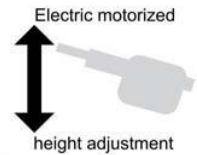

Model 2895-00XLE  
with options:  
Movable on twin castors with  
central lock (046), paper roll  
holder (010F/K)



### Basic equipment:

- Lying surface 4-sectioned
- Multifunctional and independent adjustment of upholstery sections
- Possibilities of positioning the surface parts: Roof shape, Trendelenburg, reverse Trendelenburg
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 4-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 420 mm; with nose opening and upholstery part to close
- Middle part stepless up- (+45°) and downwards (-15°) adjustable by gas spring
- Length middle part: 550 mm
- Foot part 2-sectioned, split transversal:  
Upper foot segment stepless up- (+50°) and downwards (-15°) adjustable by gas spring  
Lower foot segment stepless up- (+35°) and downwards (-30°) adjustable by gas spring
- Length foot part: Upper foot segment: 550 mm | Lower foot segment: 420 mm
- Thickness of upholstery: 60 mm

### Technical data:


- Electronical height adjustment: 570 mm to 1010 mm in 22 sec.  
by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

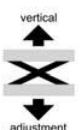

## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 077N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock



Art. no. 049  
Movable on twin castors, Ø 100 mm  
with central lock



Art. no. 044N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering



### Additional motor control options



Art. no. 090E  
Hand switch (low voltage),  
with magnet on the back to fasten



Art. no. 091E  
Foot switch (low voltage)



Art. no. 092E  
1 foot switch fixation mounted at  
underframe

### Additional safety options



Art. no. B1997E  
Foot switch cover,  
made of brushed VA-steel

### Paper roll holder



Art. no. 010F/K  
Paper roll holder, mounted under head  
part (not protruding), for width 700 mm

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery  
part to close

### Upholstery options

Art. no. B1770  
Soft covering for pleasant resting comfort,  
for width 700 mm, thickness 80 mm

Art. no. B1772  
Upholstery medium-soft (cold foam cover),  
for width 700 mm, thickness 70 mm


# Series 2800XL

## 4-sectioned therapy table with hydraulic height adjustment

**Model 2895-00XL/H white powder coated (RAL 9010)**

**Model 2895-03XL/H white aluminium powder coated (RAL 9006)**

**Model 2895-04XL/H grey aluminium powder coated (RAL 9007)**

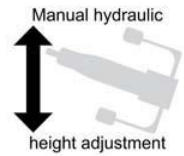

Under part in case of hydraulic height adjustment,  
with options:  
Movable on twin castors with central lock (046)

without  
picture

Model 2895-00XLE (electric motorized height adjustment) with options:  
Movable on twin castors with central lock (046),  
paper roll holder (010F/K)



### Basic equipment:

- Lying surface 4-sectioned
- Multifunctional and independent adjustment of upholstery sections
- Possibilities of positioning the surface parts: Roof shape, Trendelenburg, reverse Trendelenburg
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 4-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 420 mm; with nose opening and upholstery part to close
- Middle part stepless up- (+45°) and downwards (-15°) adjustable by gas spring
- Length middle part: 550 mm
- Foot part 2-sectioned, split transversal:
  - Upper foot segment stepless up- (+50°) and downwards (-15°) adjustable by gas spring
  - Lower foot segment stepless up- (+35°) and downwards (-30°) adjustable by gas spring
- Length foot part: Upper foot segment: 550 mm | Lower foot segment: 420 mm
- Thickness of upholstery: 60 mm



### Technical data:


- Manual foot operated height adjustment: 600 mm to 1020 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

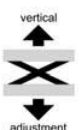



## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 077N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock



Art. no. 049  
Movable on twin castors, Ø 100 mm  
with central lock



Art. no. 044N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering



### Paper roll holder



Art. no. 010F/K  
Paper roll holder, mounted under head  
part (not protruding), for width 700 mm

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery  
part to close

### Upholstery options

Art. no. B1770  
Soft covering for pleasant resting comfort,  
for width 700 mm, thickness 80 mm

Art. no. B1772  
Upholstery medium-soft (cold foam cover),  
for width 700 mm, thickness 70 mm



# Series 2800XL

## 4-sectioned therapy table with electric motorized height adjustment

**Model 2850-00XLE** white powder coated (RAL 9010)

**Model 2850-03XLE** white aluminium powder coated (RAL 9006)

**Model 2850-04XLE** grey aluminium powder coated (RAL 9007)

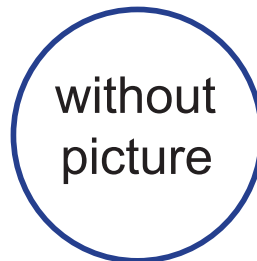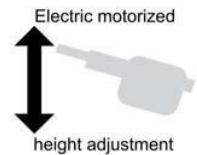

### Basic equipment:

- Lying surface 4-sectioned
- Middle part adjustable to roof position, Trendelenburg position possible
- Foot part 2-sectioned, split transversal
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 4-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part 2-sectioned, split transversal
- Complete foot section stepless adjustable (+60°) by a self-fixing system of 2 guide rods, without clearance
- Lower foot section additionally stepless adjustable (-65°) by gas spring (roof position)
- Length foot part: Upper foot section: 405 mm | Lower foot section: 485 mm
- Thickness of upholstery: 60 mm



### Technical data:


- Electronical height adjustment: 510 mm to 950 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

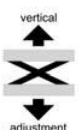

## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 040N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock,  
basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm  
with central lock, basic height +20 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering



### Additional motor control options



Art. no. 090E  
Hand switch (low voltage),  
with magnet on the back to fasten



Art. no. 091E  
Foot switch (low voltage)



Art. no. 092E  
1 foot switch fixation mounted at  
underframe



Art. no. 094E  
Foot switch rail, lateral, on both sides,  
security adjustment



Art. no. 095E  
Foot switch rail, surround,  
security adjustment

### Additional safety options



Art. no. B1997E  
Foot switch cover,  
made of brushed VA-steel

### Paper roll holder



Art. no. 010K  
Paper roll holder, mounted under head  
part (not protruding), for width 700 mm,  
not in case of 3-sectioned head part  
(Art. no. 056)

### Head part options



Art. no. 056  
Head part 3-sectioned, arm rests  
downwards adjustable from -25° to  
+45° and removable  
(art. no. 010K not possible)  
measures:  
Head part: L. 540 mm x W. 370 mm  
Arm rests: L. 500 mm x W. 160 mm

### Upholstery options

Art. no. 008  
Width of upholstery 800 mm

Art. no. B1772  
Upholstery medium-soft (cold foam  
cover), for width 700 mm,  
thickness 70 mm

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery  
part to close


# Series 2800XL

## 4-sectioned therapy table with hydraulic height adjustment

**Model 2850-00XL/H white powder coated (RAL 9010)**

**Model 2850-03XL/H white aluminium powder coated (RAL 9006)**

**Model 2850-04XL/H grey aluminium powder coated (RAL 9007)**

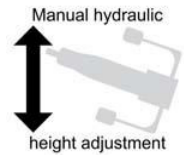

without  
picture



### Basic equipment:

- Lying surface 4-sectioned
- Middle part adjustable to roof position, Trendelenburg position possible
- Foot part 2-sectioned, split transversal
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 4-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part 2-sectioned, split transversal
- Complete foot section stepless adjustable (+60°) by a self-fixing system of 2 guide rods, without clearance
- Lower foot section additionally stepless adjustable (-65°) by gas spring (roof position)
- Length foot part: Upper foot section: 405 mm | Lower foot section: 485 mm
- Thickness of upholstery: 60 mm



### Technical data:


- Manual foot operated height adjustment: 530 mm to 950 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

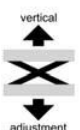



## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 040N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock,  
basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm with  
central lock, basic height +20 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering



### Paper roll holder



Art. no. 010K  
Paper roll holder, mounted under head  
part (not protruding), for width 700 mm,  
not in case of 3-sectioned head part  
(Art. no. 056)

### Head part options



Art. no. 056  
Head part 3-sectioned, arm rests down-  
wards adjustable from -25° to +45° and  
removable  
(art. no. 010K not possible)  
measures:  
Head part: L. 540 mm x W. 370 mm  
Arm rests: L. 500 mm x W. 160 mm

### Upholstery options

Art. no. 008  
Width of upholstery 800 mm

Art. no. B1772  
Upholstery medium-soft (cold foam  
cover), for width 700 mm,  
thickness 70 mm

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery  
part to close


# Series 2800XL

## 4-sectioned therapy table with electric motorized height adjustment

Model 2860-00XLE white powder coated (RAL 9010)

Model 2860-03XLE white aluminium powder coated (RAL 9006)

Model 2860-04XLE grey aluminium powder coated (RAL 9007)

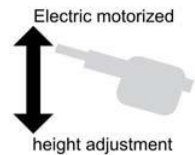

All pictures:  
Model 2860-00XLE  
in basic equipment



### Basic equipment:

- Lying surface 4-sectioned
- Middle part adjustable to roof position, Trendelenburg position possible
- Foot part 2-sectioned, split longitudinal
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 4-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part 2-sectioned, split longitudinal | Length foot part: 940 mm
- Every foot section stepless adjustable (+65°) by a self-fixing system of 2 guide rods, without clearance (foot sections independently adjustable)
- Thickness of upholstery: 60 mm



### Technical data:


- Electronical height adjustment: 510 mm to 950 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

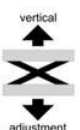

## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 040N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock,  
basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm  
with central lock, basic height +20 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering



### Additional motor control options



Art. no. 090E  
Hand switch (low voltage),  
with magnet on the back to fasten



Art. no. 091E  
Foot switch (low voltage)



Art. no. 092E  
1 foot switch fixation mounted at  
underframe



Art. no. 094E  
Foot switch rail, lateral, on both sides,  
security adjustment



Art. no. 095E  
Foot switch rail, surround,  
security adjustment

### Additional safety options



Art. no. B1997E  
Foot switch cover,  
made of brushed VA-steel

### Paper roll holder



Art. no. 010K  
Paper roll holder, mounted under head  
part (not protruding), for width 700 mm,  
not in case of 3-sectioned head part  
(Art. no. 056)

### Head part options



Art. no. 056  
Head part 3-sectioned, arm rests  
downwards adjustable from -25° to  
+45° and removable  
(art. no. 010K not possible)  
measures:  
Head part: L. 540 mm x W. 370 mm  
Arm rests: L. 500 mm x W. 160 mm

### Upholstery options

Art. no. 008  
Width of upholstery 800 mm

Art. no. B1772  
Upholstery medium-soft (cold foam  
cover), for width 700 mm,  
thickness 70 mm

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery  
part to close


## Series 2800XL

### 4-sectioned therapy table with hydraulic height adjustment

**Model 2860-00XL/H white powder coated (RAL 9010)**

**Model 2860-03XL/H white aluminium powder coated (RAL 9006)**

**Model 2860-04XL/H grey aluminium powder coated (RAL 9007)**

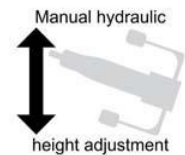

without  
picture



#### Basic equipment:

- Lying surface 4-sectioned
- Middle part adjustable to roof position, Trendelenburg position possible
- Foot part 2-sectioned, split longitudinal
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



#### Upholstery:

Lying surface 4-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part 2-sectioned, split longitudinal | Length foot part: 940 mm
- Every foot section stepless adjustable (+65°) by a self-fixing system of 2 guide rods, without clearance (foot sections independently adjustable)
- Thickness of upholstery: 60 mm



#### Technical data:


- Manual foot operated height adjustment: 530 mm to 950 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



#### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

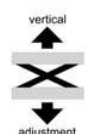



Available as an option (see pricelist for additional options)

### Mobility




Art. no. 040N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock,  
basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm with  
central lock, basic height +20 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering

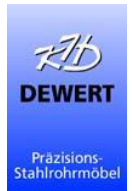

### Paper roll holder



Art. no. 010K  
Paper roll holder, mounted under head  
part (not protruding), for width 700 mm,  
not in case of 3-sectioned head part  
(Art. no. 056)

### Head part options



Art. no. 056  
Head part 3-sectioned, arm rests down-  
wards adjustable from -25° to +45° and  
removable  
(art. no. 010K not possible)  
measures:  
Head part: L. 540 mm x W. 370 mm  
Arm rests: L. 500 mm x W. 160 mm

### Upholstery options

Art. no. 008  
Width of upholstery 800 mm

Art. no. B1772  
Upholstery medium-soft (cold foam  
cover), for width 700 mm,  
thickness 70 mm

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery  
part to close


# Series 2800XL

## 6-sectioned therapy table with electric motorized height adjustment

**Model 2880-00XLE** white powder coated (RAL 9010)

**Model 2880-03XLE** white aluminium powder coated (RAL 9006)

**Model 2880-04XLE** grey aluminium powder coated (RAL 9007)

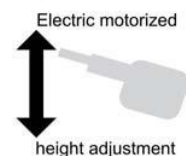

All pictures:  
Model 2880-00E  
with options:  
Wheel-lifting/fixing mechanism (045N), surrounding foot control (095E)



### Basic equipment:

- Lying surface 6-sectioned
- Middle part adjustable to roof position, Trendelenburg position possible
- Foot part 4-sectioned, split transversal and longitudinal
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 6-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part 4-sectioned, split transversal and longitudinal
- Every foot section stepless adjustable (+60°) by a self-fixing system of 2 guide rods, without clearance (foot sections independently adjustable)
- Lower foot sections additionally stepless adjustable (-65°) by gas spring
- Length foot part: Upper foot section: 405 mm | Lower foot section: 485 mm
- Thickness of upholstery: 60 mm



### Technical data:


- Electronical height adjustment: 510 mm to 950 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

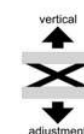

## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 040N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock,  
basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm with  
central lock, basic height +20 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering



### Additional motor control options



Art. no. 090E  
Hand switch (low voltage),  
with magnet on the back to fasten



Art. no. 091E  
Foot switch (low voltage)



Art. no. 092E  
1 foot switch fixation mounted at  
underframe



Art. no. 094E  
Foot switch rail, lateral, on both sides,  
security adjustment



Art. no. 095E  
Foot switch rail, surround,  
security adjustment



Art. no. 065E  
Electric motorized adjustment to roof  
position

### Additional safety options



Art. no. B1997E  
Foot switch cover,  
made of brushed VA-steel

### Paper roll holder



Art. no. 010K  
Paper roll holder, mounted under head  
part (not protruding), for width 700 mm,  
not in case of 3-sectioned head part  
(Art. no. 056)

### Nose opening




Art. no. 1615-EF  
Loss of nose opening and upholstery  
part to close

### Head part options



Art. no. 056  
Head part 3-sectioned, arm rests down-  
wards adjustable from -25° to +45° and  
removable  
(art. no. 010K not possible)  
measures:  
Head part: L. 540 mm x W. 370 mm  
Arm rests: L. 500 mm x W. 160 mm

### Arm rests

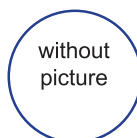

Art. no. 088  
1 arm rest, fixed, swivelling and  
removeable, additional revolving and  
tiltable

### Upholstery options

Art. no. 009  
Width of upholstery 800 mm

Art. no. B1772  
Upholstery medium-soft (cold foam  
cover), for width 700 mm,  
thickness 70 mm



# Series 2800XL

## 6-sectioned therapy table with hydraulic height adjustment

**Model 2880-00XL/H white powder coated (RAL 9010)**

**Model 2880-03XL/H white aluminium powder coated (RAL 9006)**

**Model 2880-04XL/H grey aluminium powder coated (RAL 9007)**

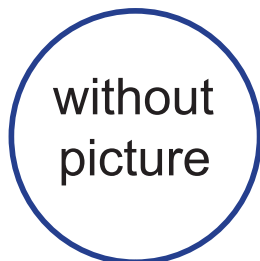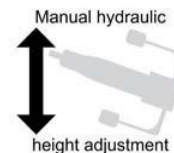

Model 2880-00E (electric motorized height adjustment) with options:  
Wheel-lifting/fixing mechanism (045N), surrounding foot control (095E)



### Basic equipment:

- Lying surface 6-sectioned
- Middle part adjustable to roof position, Trendelenburg position possible
- Foot part 4-sectioned, split transversal and longitudinal
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 6-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part 4-sectioned, split transversal and longitudinal
- Every foot section stepless adjustable (+60°) by a self-fixing system of 2 guide rods, without clearance (foot sections independently adjustable)
- Lower foot sections additionally stepless adjustable (-65°) by gas spring
- Length foot part: Upper foot section: 405 mm | Lower foot section: 485 mm
- Thickness of upholstery: 60 mm



### Technical data:


- Manual foot operated height adjustment: 590 mm to 1010 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

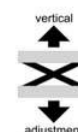



## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 040N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock,  
basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm  
with central lock, basic height +20 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering



### Paper roll holder



Art. no. 010K  
Paper roll holder, mounted under head  
part (not protruding), for width 700 mm,  
not in case of 3-sectioned head part  
(Art. no. 056)

### Head part options




Art. no. 056  
Head part 3-sectioned, arm rests  
downwards adjustable from -25° to  
+45° and removable  
(art. no. 010K not possible)  
measures:  
Head part: L. 540 mm x W. 370 mm  
Arm rests: L. 500 mm x W. 160 mm

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery  
part to close

### Arm rests

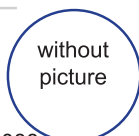

Art. no. B088  
1 arm rest, fixed, swivelling and  
removeable, additional revolving and  
tiltable

### Upholstery options

Art. no. 009  
Width of upholstery 800 mm

Art. no. B1772  
Upholstery medium-soft (cold foam  
cover), for width 700 mm,  
thickness 70 mm


# Series 4800E

## 3-sectioned therapy table with electric motorized height adjustment

**Model 4815-00E white powder coated (RAL 9010)**

**Model 4815-03E white aluminium powder coated (RAL 9006)**

**Model 4815-04E grey aluminium powder coated (RAL 9007)**

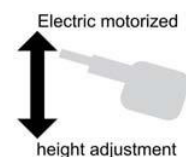

All pictures:  
Model 4815-00E  
with options:  
Movable on twin castors with  
central lock (046), surrounding  
foot control (095E), paper roll  
holder (010F)



Table in highest position



Table in lowest position



### Basic equipment:

- Lying surface 3-sectioned
- Foot part stepless adjustable (+75°) to sitting position
- Table not movable, including height adjustment screws to compensate for uneven ground
- High adjustment range
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+45°) and downwards (-35°) adjustable by gas spring
- Length head part: 500 mm; with nose opening and upholstery part to close
- Middle part: fixed
- Length middle part: 530 mm
- Foot part stepless upwards adjustable (+75°) to sitting position by gas spring
- Length foot part: 940 mm
- Thickness of upholstery: 60 mm



### Technical data:


- Electronical height adjustment: 450 mm to 1000 mm in 18 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1970 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

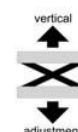

## Available as an option (see pricelist for additional options)

### Mobility




Art. no. 040N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock,  
basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm  
with central lock, basic height +20 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering

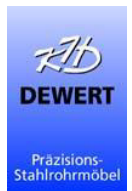

### Additional motor control options



Art. no. 090E  
Hand switch (low voltage),  
with magnet on the back to fasten



Art. no. 091E  
Foot switch (low voltage)



Art. no. 092E  
1 foot switch fixation mounted at  
underframe

### Additional safety options



Art. no. 094E  
Foot switch rail, lateral, on both sides,  
security adjustment



Art. no. 095E  
Foot switch rail, surround,  
security adjustment



Art. no. B1997E  
Foot switch cover,  
made of brushed VA-steel

### Paper roll holder



Art. no. 010K  
Paper roll holder, mounted under head  
part (not protruding), for width 700 mm,  
not in case of 3-sectioned head part  
(Art. no. 056)




Art. no. 010F  
Paper roll holder, mounted at foot part  
(protruding: 70 mm), for width 700 mm



Art. no. 055  
Head part 3-sectioned, arm rests  
downwards adjustable from -25° to  
+45° and removable  
(art. no. 010K not possible)  
measures:  
Head part: L. 500 mm x W. 370 mm  
Arm rests: L. 420 mm x W. 160 mm

### Arm rests

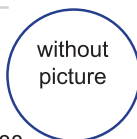

Art. no. 088  
1 arm rest, fixed, swivelling and  
removeable, additional revolving and  
tiltable

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery  
part to close

### Upholstery options

Art. no. 007  
Width of upholstery 800 mm

Art. no. B1770  
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. B1772  
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm


# Series 4800E

## 3-sectioned therapy table with electric motorized height adjustment

**Model 4840-00E white powder coated (RAL 9010)**

**Model 4840-03E white aluminium powder coated (RAL 9006)**

**Model 4840-04E grey aluminium powder coated (RAL 9007)**

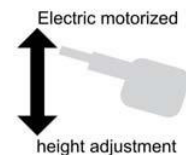

All pictures:  
Model 4840-00E  
with options:  
Wheel-lifting/fixing mechanism  
(045N),  
surrounding foot control (095E),  
head part 3-sectioned (055-10)



### Basic equipment:

- Lying surface 3-sectioned
- Multifunctional and independent adjustment of upholstery sections
- Possibilities of positioning the surface parts: Roof shape, Trendelenburg, reverse Trendelenburg
- Table not movable, including height adjustment screws to compensate for uneven ground
- High adjustment range
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+ 45°) and downwards (-35°) adjustable by gas spring
- Length head part: 450 mm
- Head part with nose opening and upholstery part to close
- Middle part adjustable to roof position by gas spring
- Length middlepart: 530 mm
- Foot part stepless adjustable (+65°) to sitting position by 2 gas springs
- Length foot part: 1020 mm
- Thickness of upholstery: 60 mm



### Technical data:

- Electronical height adjustment: 510 mm to 1030 mm in 18 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 2000 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology





## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 040N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock,  
basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm  
with central lock, basic height +20 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering



### Additional motor control options



Art. no. 090E  
Hand switch (low voltage),  
with magnet on the back to fasten



Art. no. 091E  
Foot switch (low voltage)



Art. no. 092E  
1 foot switch fixation mounted at  
underframe

### Additional safety options



Art. no. 094E  
Foot switch rail, lateral, on both sides,  
security adjustment



Art. no. 095E  
Foot switch rail, surround,  
security adjustment



Art. no. B1997E  
Foot switch cover,  
made of brushed VA-steel

### Paper roll holder



Art. no. 010K  
Paper roll holder, mounted under head  
part (not protruding), for width 700 mm,  
not in case of 3-sectioned head part  
(Art. no. 056)




Art. no. 010F  
Paper roll holder, mounted at foot part  
(protruding: 70 mm), for width 700 mm



Art. no. 055-10  
Head part 3-sectioned, arm rests  
downwards adjustable from -25° to  
+45° and removable  
(art. no. 010K not possible)  
measures:  
Head part: L. 500 mm x W. 370 mm  
Arm rests: L. 420 mm x W. 160 mm

### Arm rests

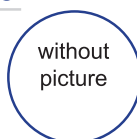

Art. no. B088  
1 arm rest, fixed, swivelling and  
removeable, additional revolving and  
tiltable

### Nose opening



Art. no. 1615-EF  
Loss of nose opening and upholstery  
part to close

### Upholstery options

Art. no. 007  
Width of upholstery 800 mm

Art. no. B1772  
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

# Series 2800XL


**Gymnastic table for therapy acc. to Bobath and Vojta  
with electric motorized height adjustment**

**Model 2800-00XLE\* white powder coated (RAL 9010)**

**Model 2800-03XLE\* white aluminium powder coated (RAL 9006)**

**Model 2800-04XLE\* grey aluminium powder coated (RAL 9007)**

\* Exact name of model depending on width of lying surface (see technical data)

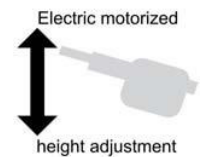

Model 2804-00XLE  
with options:  
Wheel-lifting/fixing mechanism  
(045N), surrounding foot control  
(095E), head part stepless  
adjustable by gas spring (063)



## Basic equipment:

- Lying surface 1-sectioned
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Table in lowest position



Both pictures.:  
Model 2800-00XLE with  
option:  
Movable on twin  
castors with central lock  
(046)

## Upholstery:

Lying surface 1-sectioned:

- Length: 1950 mm
- Optional: Head part
- Thickness of upholstery: 60 mm



## Technical data:


- Electronical height adjustment: 480 mm to 920 mm in 22 sec.  
by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: **Model 2800XLE: 1000 mm** | **Model 2802XLE: 1100 mm** | **Model 2804XLE: 1200 mm**
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



## Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

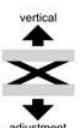

## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 040N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock,  
basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm  
with central lock, basic height +20 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering



### Additional motor control options



Art. no. 090E  
Hand switch (low voltage),  
with magnet on the back to fasten



Art. no. 091E  
Foot switch (low voltage)



Art. no. 092E  
1 foot switch fixation mounted at  
underframe



Art. no. 094E  
Foot switch rail, lateral, on both sides,  
security adjustment



Art. no. 095E  
Foot switch rail, surround,  
security adjustment

### Additional safety options



Art. no. B1997E  
Foot switch cover,  
made of brushed VA-steel

### Paper roll holder



Art. no. 016K  
Paper roll holder, mounted under head  
part (not protruding), for width 1000 mm

### Head part options



Art. no. 063  
Head part, upwards and downwards  
adjustable by gas spring

# Series 4800E


## Gymnastic table for therapy acc. to Bobath and Vojta with electric motorized height adjustment

**Model 4800-00E\*** white powder coated (RAL 9010)

**Model 4800-03E\*** white aluminium powder coated (RAL 9006)

**Model 4800-04E\*** grey aluminium powder coated (RAL 9007)

\* Exact name of model depending on width of lying surface (see technical data)

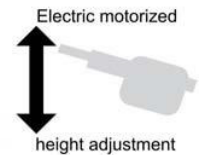

All pictures:  
Model 4804-00E  
with option:  
Wheel-lifting/fixing mechanism  
(045N)



### Basic equipment:

- Lying surface 1-sectioned
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- High adjustment range
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Table in lowest position



Table in highest position

### Upholstery:

Lying surface 1-sectioned:

- Length: 1950 mm
- Optional: Head part
- Thickness of upholstery: 60 mm



### Technical data:


- Electronical height adjustment: 420 mm to 970 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: **Model 4800E:** 1000 mm | **Model 4802E:** 1100 mm | **Model 4804E:** 1200 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

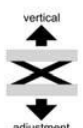



## Available as an option (see pricelist for additional options)

### Mobility



Art. no. 040N  
Movable on comfort twin castors  
Ø 100 mm, each with total lock,  
basic height -15 mm



Art. no. 046  
Movable on twin castors, Ø 100 mm  
with central lock, basic height +20 mm



Art. no. 045N  
Wheel-lifting/fixing mechanism  
with 4 foot levers for controlled  
lowering



### Additional motor control options



Art. no. 090E  
Hand switch (low voltage),  
with magnet on the back to fasten



Art. no. 091E  
Foot switch (low voltage)



Art. no. 092E  
1 foot switch fixation mounted at  
underframe



Art. no. 094E  
Foot switch rail, lateral, on both sides,  
security adjustment



Art. no. 095E  
Foot switch rail, surround,  
security adjustment

### Additional safety options



Art. no. B1997E  
Foot switch cover,  
made of brushed VA-steel

### Paper roll holder



Art. no. 016K  
Paper roll holder, mounted under head  
part (not protruding), for width 1000 mm

### Head part options



Art. no. 063  
Head part, upwards and downwards  
adjustable by gas spring

K.H. Dewert GmbH  
Vollmestraße 7  
33649 Bielefeld  
Germany

Post office 14 40 28  
33634 Bielefeld  
Germany

Phone + 49 (0) 5 21 / 4 00 27-0  
Fax + 49 (0) 5 21 / 4 00 27-27

Internet: [www.khdewert.de](http://www.khdewert.de)  
E-Mail: [info@khdewert.de](mailto:info@khdewert.de)  
YouTube: K.H. Dewert



05.2021\_01

## THERAPY TABLES

### + Additional equipment

#### Cushions



Art. no. 1211-..\*  
Wedge shaped bolster,  
measures (L. x W. x H.):  
400 x 400 x 100 - 10 mm

\*Add colour number



Art. no. 1113-..\*  
Half roll bolster,  
measures (L. x W. x H.):  
570 x 200 x 95 mm



Art. no. 1111-..\*  
Neck roll,  
measures:  
L .570 x Ø 150 mm



Art. no. 1112-..\*  
Knee roll,  
measures:  
L. 570 x Ø 200 mm



Art. no. 1314-..\*  
Cushion, velcro fixation,  
ticking with  
polyester filling,  
measures (L. x W.):  
400 x 400 mm

#### Upholstery protection mat



Art. no. 1890  
Protection mat, made of PVC,  
thickness: 0.8 mm,  
measures: 700 x 600 mm

