


MULTIPURPOSE TABLES
FOR EXAMINATION, REKTOSCOPY UND GYNAECOLOGY

SERIES **2000XL**

2021.1

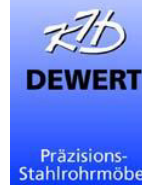

... WITH 3-SECTIONED LYING SURFACE



... WITH 4-SECTIONED LYING SURFACE



... WITH INCLINATION ADJUSTMENT



TABLES FOR SURGERY AND HOSPITAL

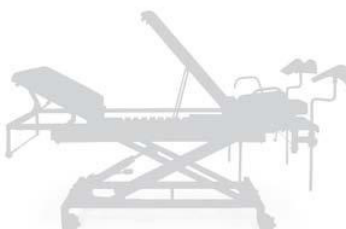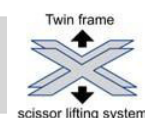

Series 2000XL - Tables for gynaecology and rectoscopy

Good reasons for a multipurpose table of series 2000XL

Content

Good reasons for a multipurpose table of series 2000XL 3

Multipurpose tables for examination, rectoscopy and gynaecology

4-sectioned multipurpose tables:

- Model 2640XLE 6
- Model 2640XL/H 8

3-sectioned multipurpose tables:

- Model 2645XLE 10
- Model 2645XL/H 12

3-sectioned multipurpose tables with inclination adjustment:

- Model 2648XLE 14





Series 2000XL - Tables for gynaecology and rectoscopy

Good reasons for a multipurpose table of series 2000XL




Electric motorized height adjustment

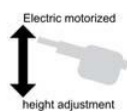

- Electronical height adjustment by AC-motor with a pressure capacity of 10000 N
- Constant adjustment speed even in case of high load
- EPR-safety power cable with strain relief, length: 4 m
- **Integrated safety devices (according to the newest advice of BfArM)** 
 - Automatic blocking of electromotive drive
 - Active freewheeling downwards, automatic declutch in case of resistance
- Load capacity: 220 kg



Hydraulic height adjustment

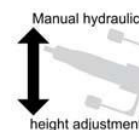

- Manual foot operated height adjustment by a hydraulic pump with a pressure capacity of 10000 N
- Pedals for foot-controlled height positioning on both sides of the table
- Load capacity: 220 kg



Dewert-hygiene-equipment Enables optimal cleaning and disinfection

- Artificial leather underneath the upholstery
- Vent holes underneath the upholstery, closed by special caps with valve function
- Easy hygiene because of an open construction: good accessibility of components, no plastic coverings which are breakable in the ageing process
- Option: Innovative artificial leather Skai® Toronto EN with staynu



Series 2000XL - Tables for gynaecology and rectoscopy

Good reasons for a multipurpose table of series 2000XL




Durable ✓




Robust ✓



Easy-care ✓

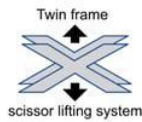

**Twin frame
scissor lifting system** ✓

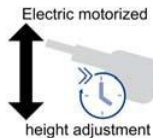

**High
adjustment speed** ✓



**Proven
Dewert-quality** ✓



Technology:



- Extremely stable twin frame scissor lifting system, with high torsional stiffness due to thick-walled square tubes with big dimensions, also for under and upper frame and scissor


- High stability due to massive 20 mm axes



- Widely-spaced feet and a high deadweight guarantee stability
- Robust frame construction which is proven in demanding everyday life of hospitals for years



- Abrasion-resistant guide system at under and upper frame due to the use of roller bearings

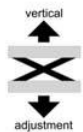


Upholstery:

- Thickness of upholstery: 60 mm
- Multiple glued, high-strength laminated wood panel (thickness: 18 mm) no breaking chipboards in case of pressure on the edges
- Two-layer homogenous polyurethane foam, density 46 and 75, Ökotex-certified
- 2 different variants for artificial leather eligible for upholstery covering:
 - Skai® Pandoria Plus with velvet-smooth surface
 - Skai® Toronto EN with innovative staynu-technology



Series 2000XL - Tables for gynaecology and rectoscopy

Good reasons for a multipurpose table of series 2000XL

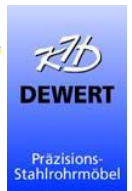

... with 3-sectioned or 4-sectioned lying surface:

- Head part 30° positive adjustable by 2 metal rachets
- Middle part 50° positive adjustable by rack system
- Seat section with half-round cut-out, manually inclinable to +20°
- 4-sectioned lying surface additionally:
Rectoscopy insert to lift up, width: 100 mm



- 3-sectioned lying surface and inclination adjustment additionally:
Electric motorized inclination adjustment to -15°, operation by hand switch



Equipment:

- Hand grips on both sides of the seat section
- Drawable basin (stainless steel) under seat basin, support: HPL panel
- Movable on 4 comfort twin castors, Ø 100 mm, each with lock
- Clamps for leg holders, variants of leg supports available
- Dewert-hygiene-equipment




You will find more options on the following pages.

Series 2000XL

4-sectioned multipurpose table with electric motorized height adjustment

... for examination, rectoscopy, gynaecology ...

Model 2640-01XLE top frame chromium-painted, underpart white powder coated (RAL 9010)

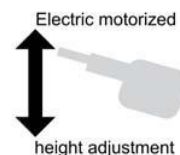

All pictures:
Modell 2640-01XLE
with option:
Paper roll holder (1265K);
big picture additionally: leg
supports acc. to Goepel
(1520)



Basic equipment:

- Movable on comfort twin castors, Ø 100 mm, each with lock
- Hand grips on both sides of the seat section
- Drawable basin (stainless steel) under seat, basin support: HPL panel
- Including clamps for leg holders (leg holders = available as an option)
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 4-sectioned:

- Head part 30° positive adjustable by 2 metal rachets
- Length head part: 550 mm
- Middle part 50° positive adjustable by rack system
- Length middle part: 930 mm
- Rectoscopy insert to lift up and fixable, width: 100 mm
- Seat section with half-round cut-out, manually inclinable (20°)
- Length seat section: 400 mm
- Thickness of upholstery: 80 mm



Technical data:

- Electronical height adjustment: 620 mm to 1060 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 650 mm
- Length of lying surface: 1980 mm
- Load capacity: 220 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1182
Without basic mobility



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close



Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe

Additional safety options



Art. no. 1950E
Foot switch rail, lateral, on both sides, security adjustment



Art. no. B1997E
Foot switch cover, made of brushed VA-steel

Paper roll holder



Art. no. 1265K
Paper roll holder, mounted under head part (not protruding), for with of lying surface 700 mm

Side guards



Art. no. 1040
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm

Head part options



Art. no. 1655
Head part stepless adjustable by gas spring from 0° to +60°

Add-on elements



Art. no. 1410
1 piece of side rail, length: 350 mm



Art. no. 1430
1 attachment clamp for side rail, suitable for pipe Ø 18 mm



Art. no. 1450
1 Infusion pole for side rail, without attachment clamp

Leg supports



Art. no. 1510
1 pair of leg supports



Art. no. 1520
1 pair of leg supports according to Goepel

Upholstery options

Art. no. 1770
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm


Art. no. 1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

Series 2000XL

4-sectioned multipurpose table with hydraulic height adjustment

... for examination, rectoscopy, gynaecology ...

Model 2640-01XL/H top frame chromium-painted, underpart white powder coated (RAL 9010)

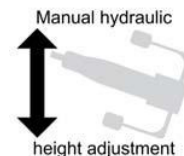

2 upper pictures:
Model 2640-00XL/H
with option:
Paper roll holder (1265K),
side rail (1410),
leg supports (1510)



Basic equipment:

- Movable on comfort twin castors, Ø 100 mm, each with lock
- Hand grips on both sides of the seat section
- Drawable basin (stainless steel) under seat, basin support: HPL panel
- Including clamps for leg holders (leg holders = available as an option)
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-Hygiene-Ausstattung



Upholstery:

Lying surface 4-sectioned:

- Head part 30° positive adjustable by 2 metal rachets
- Length head part: 550 mm
- Middle part 50° positive adjustable by rack system
- Length middle part: 930 mm
- Rectoscopy insert to lift up and fixable, width: 100 mm
- Seat section with half-round cut-out, manually inclinable (20°)
- Length seat section: 400 mm
- Thickness of upholstery: 80 mm



Model 2640-00XL/H with option:
Paper roll holder, side rail (1410)

Technical data:

- Manual foot operated height adjustment: 640 mm to 1060 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: 650 mm
- Length of lying surface: 1980 mm
- Load capacity: 220 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1182
Without basic mobility



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Paper roll holder



Art. no. 1265K
Paper roll holder, mounted under head part (not protruding), for width of lying surface 700 mm

Side guards



Art. no. 1040
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm

Head part options



Art. no. 1655
Head part stepless adjustable by gas spring from 0° to +60°

Add-on elements



Art. no. 1410
1 piece of side rail, length: 350 mm



Art. no. 1430
1 attachment clamp for side rail, suitable for pipe Ø 18 mm



Art. no. 1450
1 Infusion pole for side rail, without attachment clamp

Leg supports



Art. no. 1510
1 pair of leg supports



Art. no. 1520
1 pair of leg supports according to Goepel

Upholstery options



Art. no. 1770
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. 1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

Series 2000XL

3-sectioned multipurpose table with electric motorized height adjustment ... for examination, rectoscopy, gynaecology ...

**Model 2645-01XLE top frame chromium-painted,
underpart white powder coated (RAL 9010)**

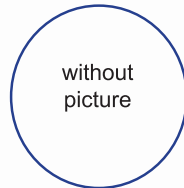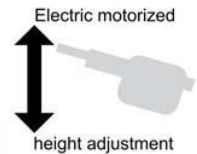

Model 2640-01XLE with option: Paper roll holder (1265K)

Model 2645-01XLE differs from Model 2640-01XLE
(see picture):

Lying surface 3-sectioned, no rectoscopy insert



Basic equipment:

- Movable on comfort twin castors, Ø 100 mm, each with lock
- Hand grips on both sides of the seat section
- Drawable basin (stainless steel) under seat, basin support: HPL panel
- Including clamps for leg holders (leg holders = available as an option)
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

lying surface 3-teilig:

- Head part 30° positive adjustable by 2 metal rachets
- Length head part: 550 mm
- Middle part 50° positive adjustable by rack system
- Length middle part: 930 mm
- Seat section with half-round cut-out, manually inclinable (20°)
- Length seat section: 450 mm
- Thickness of upholstery: 80 mm



Technical data:



- Electronical height adjustment: 620 mm to 1060 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 650 mm
- Length of lying surface: 1980 mm
- Load capacity: 220 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

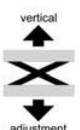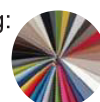

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1182
Without basic mobility



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close



Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe

Additional safety options



Art. no. 1950E
Foot switch rail, lateral, on both sides, security adjustment



Art. no. B1997E
Foot switch cover, made of brushed VA-steel

Paper roll holder



Art. no. 1265K
Paper roll holder, mounted under head part (not protruding), for width of lying surface 700 mm

Side guards



Art. no. 1040
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm

Head part options



Art. no. 1655
Head part stepless adjustable by gas spring from 0° to +60°

Add-on elements



Art. no. 1410
1 piece of side rail, length: 350 mm



Art. no. 1430
1 attachment clamp for side rail, suitable for pipe Ø 18 mm



Art. no. 1450
1 Infusion pole for side rail, without attachment clamp

Leg supports



Art. no. 1510
1 pair of leg supports



Art. no. 1520
1 pair of leg supports according to Goepel

Upholstery options

Art. no. 1770
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm


Art. no. 1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

Series 2000XL

3-sectioned multipurpose table with hydraulic height adjustment

... for examination, rectoscopy, gynaecology ...

Model 2645-01XL/H top frame chromium-painted, underpart white powder coated (RAL 9010)

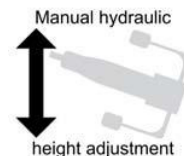

without picture



Model 2645-01XL/H differs from Model 2640-01XL/H (see picture):

Lying surface 3-sectioned, no rectoscopy insert

Both pictures:
Model 2640-00XL/H with options:
Paper roll holder (1265K), side rail (1410),
leg supports (1510)



Basic equipment:

- Movable on comfort twin castors, Ø 100 mm, each with lock
- Hand grips on both sides of the seat section
- Drawable basin (stainless steel) under seat, basin support: HPL panel
- Including clamps for leg holders (leg holders = available as an option)
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part 30° positive adjustable by 2 metal rachets
- Length head part: 550 mm
- Middle part 50° positive adjustable by rack system
- Length middle part: 930 mm
- Seat section with half-round cutout, manually inclinable (20°)
- Length seat section: 450 mm
- Thickness of upholstery: 80 mm

Technical data:

- Manual foot operated height adjustment: 640 mm to 1060 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: 650 mm
- Length of lying surface: 1980 mm
- Load capacity: 220 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1182
Without basic mobility



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Paper roll holder



Art. no. 1265K
Paper roll holder, mounted under head part (not protruding), for width of lying surface 700 mm

Side guards



Art. no. 1040
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm

Head part options



Art. no. 1655
Head part stepless adjustable by gas spring from 0° to +60°

Add-on elements



Art. no. 1410
1 piece of side rail, length: 350 mm



Art. no. 1430
1 attachment clamp for side rail, suitable for pipe Ø 18 mm



Art. no. 1450
1 Infusion pole for side rail, without attachment clamp

Leg supports



Art. no. 1510
1 pair of leg supports



Art. no. 1520
1 pair of leg supports according to Goepel

Upholstery options

Art. no. 1770
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. 1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

Series 2000XL

3-sectioned multipurpose table with electric motorized height and inclination adjustment ... for examination, rectoscopy, gynaecology ...

Model 2648-01XLE top frame chromium-painted, underpart white powder coated (RAL 9010)

Electric motorized
height adjustment



Model 2648-01XLE
with options:
Loss of basic mobility (1182);
upper and lower pictures
additionally: Leg supports (1510)



Basic equipment:

- Movable on comfort twin castors, Ø 100 mm, each with lock
- Hand grips on both sides of the seat section
- Drawable basin (stainless steel) under seat, basin support: HPL panel
- Including clamps for leg holders (leg holders = available as an option)
- Electromotive drive control of height adjustment by foot switch
- Electromotive drive control of inclination adjustment by hand switch, inclination adjustment to -15°
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part 30° positive adjustable by 2 metal rachets
- Length head part: 550 mm
- Middle part 50° positive adjustable by rack system
- Length middle part: 930 mm
- Seat section with half-round cutout, manually inclinable (20°)
- Length seat section: 450 mm
- Thickness of upholstery: 80 mm



Technical data:

- Electronical height adjustment: 630 mm to 1070 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 650 mm
- Length of lying surface: 1950 mm
- Load capacity: 220 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1182
Without basic mobility



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close



Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe



Art. no. 1950E
Foot switch rail, lateral, on both sides, security adjustment

Additional safety options



Art. no. B1997E
Foot switch cover, made of brushed VA-steel

Paper roll holder



Art. no. 1265K
Paper roll holder, mounted under head part (not protruding), for with of lying surface 700 mm

Side guards



Art. no. 1040
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm

Head part options



Art. no. 1655
Head part stepless adjustable by gas spring from 0° to +60°

Add-on elements



Art. no. 1410
1 piece of side rail, length: 350 mm



Art. no. 1430
1 attachment clamp for side rail, suitable for pipe Ø 18 mm



Art. no. 1450
1 Infusion pole for side rail, without attachment clamp

Leg supports



Art. no. 1510
1 pair of leg supports



Art. no. 1520
1 pair of leg supports according to Goepel

Upholstery options

Art. no. 1770
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. 1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm



K.H. Dewert GmbH
Vollmestraße 7
33649 Bielefeld
Germany

Post office 14 40 28
33634 Bielefeld
Germany

Phone + 49 (0) 5 21 / 4 00 27-0
Fax + 49 (0) 5 21 / 4 00 27-27

Internet: www.khdewert.de
E-Mail: info@khdewert.de
YouTube: K.H. Dewert

MULTIPRPOSE TABLES FOR EXAMINATION, RECTOSCOPY, GYNAECOLOGY

05.2021_02

TABLES FOR SURGERY AND HOSPITAL

SERIES **2000XL**

