

EXAMINATION AND THERAPY TABLES

RD
DEWERT

Präzisions-
Stahlrohrmöbel

SERIES **DX1**

2021.2



DURABLE

ROBUST

EASY-CARE



SERIES 1000



SERIES 1030



SERIES 1050



SERIES 1090

Series DX1


The models of series DX1 ensure optimal support in examination and treatment.

A proven and reliable height adjustment is either effected by a robust hydraulic pump with an optimal ratio or a high-performance freewheeling electric motor which corresponds to the newest advice of BfArM.

Equipped with our proven upholstery made of laminated wood with homogenous polyurethane foam, the models of series DX1 enable a pleasant resting comfort.

Furthermore, every table is equipped with the Dewert-hygiene-package, which guarantees optimal cleaning and disinfection of the upholstery as well as of the frame construction.

Therefore, the high quality-models of series DX1 have best prerequisites for long-lasting daily practice use.

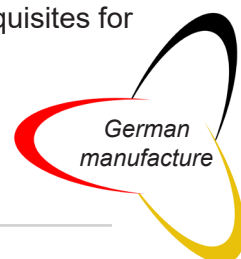

Content

Good reasons for a table of series DX1	3
--	---

Versatile tables for examination, treatment, therapy:

2-sectioned:

• Model 1000E with electric motorized height adjustment	6
• Model 1000/H with hydraulic height adjustment	8
• Model 1030E with long head part and electric motorized height adjustment	10
• Model 1030/H with long head part and hydraulic height adjustment	12

3-sectioned:

• Model 1050E with electric motorized height adjustment	14
• Model 1050/H with hydraulic height adjustment	16

Recovery and patient transport tables:

2-sectioned:

• Model 1090/H with hydraulic height adjustment	18
• Model 1093/H with long head part and hydraulic height adjustment	20

3-sectioned:

• Model 1095/H with hydraulic height adjustment	22
--	----

Additional equipment:

Cushions and upholstery protection	23
--	----




Series DX1

Good reasons for a table of series DX1




Electric motorized height adjustment

- Electronical height adjustment by AC-motor with a pressure capacity of 10000 N
- Constant adjustment speed even in case of high load
- EPR-safety power cable with strain relief, length: 4 m
- **Integrated safety devices (according to the newest advice of BfArM)**
 - Automatic blocking of electromotive drive
 - Active freewheeling downwards, automatic declutch in case of resistance
- Load capacity: 225 kg

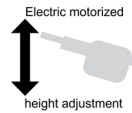

Hydraulic height adjustment

- Manual foot operated height adjustment by a hydraulic pump with a pressure capacity of 10000 N
- Pedals for foot-controlled height positioning **on both sides** of the table
- Load capacity: 225 kg

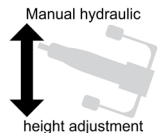

Dewert-hygiene-equipment

Enables optimal cleaning and disinfection

- Hinge lid made of identical upholstery covering
- Artificial leather underneath the upholstery
- Vent holes underneath the upholstery, closed by special caps with valve function
- Easy hygiene because of an open construction: good accessibility of components, no plastic coverings which are breakable in the ageing process
- Option: Innovative artificial leather Skai Toronto EN with staynu



Series DX1

Good reasons for a table of series DX1



Technology:



- Very stable single frame scissor lifting system, with high torsional stiffness due to thick-walled square tubes with big dimensions, also for under and upper frame and scissor
- High stability due to massive 20 mm axes



- Widely-spaced feet and a high deadweight guarantee stability
- Robust frame construction which is proven in demanding everyday life of hospitals for years



- Abrasion-resistant guide system at under and upper frame due to the use of roller bearings



Upholstery:



- Thickness of upholstery: 65 mm, recovery tables: 70 mm with a cold foam cover
- Multiple glued, high-strength laminated wood panel (thickness: 18 mm)
→ no breaking chipboards in case of pressure on the edges
- Two-layer homogenous polyurethane foam, density 46 and 75, Ökotex-certified
- 2 different variants for artificial leather eligible for upholstery covering:
 - Skai® Pandora Plus with velvet-smooth surface
 - Skai® Toronto EN with innovative staynu-technology



Series DX1

Good reasons for a table of series DX1

Lying surface 3-sectioned:



- Adjustable foot part with tubular steel frame for high torsion rigidity
- Pressure stable gas spring ensures working at the foot part without vibrations
- Operating lever to manipulate the gas spring easily accessible, positioned at the foot end
- Foot part usable as back rest for reversible positioning of the patient



Recovery and patent transport table:



- Good handling when pushing the table due to four push handles placed at the frame outside points
- Longitudinal and lateral wall protection due to four wall bumpers
- Increased lying comfort due to medium-soft upholstery with a cold foam cover
- Side guards lay down besides upholstery, below upholstery surface level, are completely made of steel tube without any breakable plastic components
- Movable on twin castors with central lock




Series DX1

2-sectioned versatile treatment tables with electric motorized height adjustment
...for examination, treatment, therapy...

Model 1000-00E white powder coated (RAL 9010)

Model 1000-03E white aluminium powder coated (RAL 9006)

Model 1000-04E silver aluminium powder coated (RAL 9007)

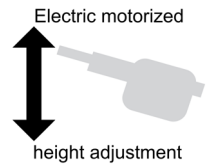

Modell 1000-00E
in basic equipment, with
hand switch



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 550 mm, integrated nose opening and upholstery part to close
- Fixed foot part
- Length foot part: 1400 mm
- Upholstery thickness: 65 mm



Modell 1000-04E mit options:
Wheel lifting/fixing mechanism (045N),
paper roll holder (010F)



Technical data:


- Electronical height adjustment: 480 - 920 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width overall: 700 mm
- Length overall: 1950 mm
- Load capacity: 225 kg, according to IEC 60601 tested with fourfold static load
Load capacity head part: 90 kg



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

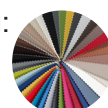

Available as an option (see pricelist for additional options)

Mobility



art. no. 040N
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



art. no. 046
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



art. no. B1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



art. no. 045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Additional motor control options



art. no. 090E
Hand switch (low voltage), with magnet on the back to fasten



art. no. 091E
Foot switch (low voltage)



art. no. 092E
1 foot switch fixation mounted at underframe



art. no. 094E
Foot switch rail, lateral, on both sides, security adjustment:
Press the switch rail down = the table is lifted
Lift the switch rail up = the table is lowered



art. no. 096E
Foot switch rail, surround, security adjustment:
Press the switch rail down = the table is lifted
Lift the switch rail up = the table is lowered



art. no. B1997E
Foot switch cover, made of brushed VA-steel

Paper roll holder



art. no. 010K/015K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm




art. no. 010F/015F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm



art. no. B1175 / B1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Head part options

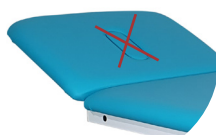

art. no. 1615-EF
Loss of nose opening and upholstery part to close



art. no. 056
Head part 3-sectioned, arm rests downwards adjustable and removable (art. no. 010K/015K not possible)

Upholstery options

art. no. 006: Width 800 mm

art. no. B1770 / B1788:
Soft covering for pleasant resting comfort, for width 700 mm/ 800 mm, thickness 80 mm

art. no. B1772 / B1782:
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Side guards



art. no. 1008
2 side guards foldable, RAL 9006 powder coated, made of steel tube, without any breakable plastic components, length: 820 mm




Series DX1

2-sectioned versatile treatment tables with hydraulic height adjustment
...for examination, treatment, therapy...

Model 1000-00/H white powder coated (RAL 9010)

Model 1000-03/H white aluminium powder coated (RAL 9006)

Model 1000-04/H silver aluminium powder coated (RAL 9007)

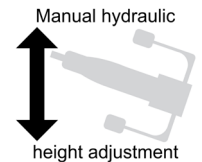

Model 1000-04/H with options:
Movable on castors with central lock (046), loss of nose opening (1615-EF)



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 550 mm, integrated nose opening and upholstery part to close
- Fixed foot part
- Length foot part: 1400 mm
- Upholstery thickness: 65 mm



Model 1000-00/H with options:
Movable on comfort twin castors, each with lock (040N), side guards (1008), loss of nose opening (1615-EF)



Technical data:

- Manual foot operated height adjustment: 470 - 920 mm by a hydraulic pump with a pressure capacity of 10000 N
 - Width overall: 700 mm
 - Length overall: 1950 mm
 - Load capacity: 225 kg, according to IEC 60601 tested with fourfold static load
- Load capacity head part: 90 kg



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with staynu-technology



Available as an option (see pricelist for additional options)

Mobility



art. no. 040N
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



art. no. 046
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



art. no. B1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



art. no. 045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Paper roll holder



art. no. 010K/015K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



art. no. 010F/015F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Side guards



art. no. 1008
2 side guards foldable, RAL 9006 powder coated, made of steel tube, without any breakable plastic components, length: 820 mm

Wall protection bumpers




art. no. B1175 / B1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Head part options



art. no. 056
Head part 3-sectioned, arm rests downwards adjustable from -25° to +45° and removable (art. no. 010K/015K not possible)
measures:
Head part: L. 550 mm x W. 380 mm
Arm rests: L. 500 mm x W. 160 mm

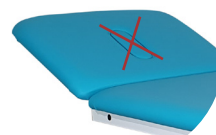

art. no. 1615-EF
Loss of nose opening and upholstery part to close

Upholstery options

art. no. 006
Width 800 mm

art. no. B1770 / B1788
Soft covering for pleasant resting comfort, for width 700 mm/ 800 mm, thickness 80 mm

art. no. B1772 / B1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm


Series DX1

2-sectioned versatile treatment tables with electric motorized height adjustment
...for examination, treatment, therapy...

Model 1030-00E white powder coated (RAL 9010)

Model 1030-03E white aluminium powder coated (RAL 9006)

Model 1030-04E silver aluminium powder coated (RAL 9007)

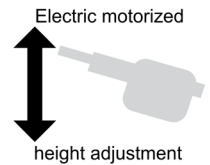

NEW



Model 1030-04E with options:
Movable on comfort twin castors,
each with lock (040N)

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Long head part with half-round design
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless adjustable from 0° to +75° by gas spring (sitting position)
- Length head part: 750 mm, without nose opening
- Fixed foot part
- Length foot part: 1200 mm
- Upholstery thickness: 65 mm



Model 1030-04E with options:
Movable on comfort twin castors, each with lock (040N)



Technical data:

- Electronical height adjustment: 520 - 960 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width overall: 700 mm
- Length overall: 1950 mm
- Load capacity: 225 kg, according to IEC 60601 tested with fourfold static load
Load capacity head part: 90 kg

kg

Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



art. no. 040N
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



art. no. 046
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



art. no. B1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



art. no. 045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Additional motor control options



art. no. 090E
Hand switch (low voltage), with magnet on the back to fasten



art. no. 091E
Foot switch (low voltage)



art. no. 092E
1 foot switch fixation mounted at underframe



art. no. 094E
Foot switch rail, lateral, on both sides, security adjustment:
Press the switch rail down = the table is lifted
Lift the switch rail up = the table is lowered



art. no. 096E
Foot switch rail, surround, security adjustment:
Press the switch rail down = the table is lifted
Lift the switch rail up = the table is lowered

Additional safety options



art. no. B1997E
Foot switch cover, made of brushed VA-steel

Paper roll holder



art. no. 010K/015K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm




art. no. 010K/F / 015K/F
Paper roll holder, mounted at foot part (not protruding), for width 700 mm / 800 mm

Wall protection bumpers



art. no. B1175 / B1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Head part option

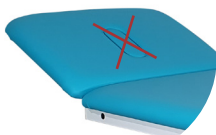

art. no. 1615-EF
Loss of nose opening and upholstery part to close

Side guards



art. no. 1008
2 side guards foldable, RAL 9006 powder coated, made of steel tube, without any breakable plastic components, length: 820 mm

Upholstery options

art. no. 006: Width 800 mm

art. no. B1770 / B1788:
Soft covering for pleasant resting comfort, for width 700 mm/ 800 mm, thickness 80 mm

art. no. B1772 / B1782:
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm


Series DX1

2-sectioned versatile treatment tables with hydraulic height adjustment
...for examination, treatment, therapy...

Model 1030-00/H white powder coated (RAL 9010)

Model 1030-03/H white aluminium powder coated (RAL 9006)

Model 1030-04/H silver aluminium powder coated (RAL 9007)

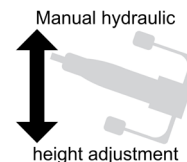

NEW



Model 1030-04/H
with options:
Movable on twin castors with
central lock (046)

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Long head part with half-round design
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless adjustable from 0° to +75° by gas spring (sitting position)
- Length head part: 750 mm, without nose opening
- Fixed foot part
- Length foot part: 1200 mm
- Upholstery thickness: 65 mm



Model 1030-04/H with options:
Movable on twin castors with central lock (046)



Technical data:


- Manual foot operated height adjustment: 510 - 960 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width overall: 700 mm
- Length overall: 1950 mm
- Load capacity: 225 kg, according to IEC 60601 tested with fourfold static load
- Load capacity head part: 90 kg



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

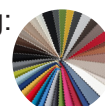

Available as an option (see pricelist for additional options)

Mobility



art. no. 040N
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



art. no. 046
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



art. no. B1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



art. no. 045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Paper roll holder



art. no. 010K/015K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



art. no. 010K/F / 015K/F
Paper roll holder, mounted at foot part (not protruding), for width 700 mm / 800 mm

Side guards



art. no. 1008
2 side guards foldable, RAL 9006 powder coated, made of steel tube, without any breakable plastic components, length: 820 mm

Wall protection bumpers



art. no. B1175 / B1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Upholstery options

art. no. 006
Width 800 mm

art. no. B1770 / B1788
Soft covering for pleasant resting comfort, for width 700 mm/ 800 mm, thickness 80 mm

art. no. B1772 / B1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm


Series DX1

3-sectioned versatile treatment tables with electric motorized height adjustment
...for examination, treatment, therapy...

Model 1050-00E white powder coated (RAL 9010)

Model 1050-03E white aluminium powder coated (RAL 9006)

Model 1050-04E silver aluminium powder coated (RAL 9007)

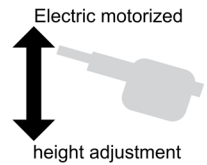

Model 1050-00E with options:
Movable on twin castors with
central lock (046), loss of nose
opening (1615-EF)



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+40°) and downwards (-30°) adjustable by gas spring
- Length head part: 550 mm, integrated nose opening and upholstery part to close
- Fixed middle part
- Length middle part: 480 mm
- Foot part stepless adjustable by gas spring (0° to +50°)
- Length foot part: 920 mm
- Upholstery thickness: 65 mm



Model 1050-00E with options:
Movable on twin castors with central lock (046), loss
of nose opening (1615-EF)



Technical data:

- Electronical height adjustment: 500 - 940 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width overall: 700 mm
- Length overall: 1950 mm
- Load capacity: 225 kg, according to IEC 60601 tested with fourfold static load
- Load capacity head part: 90 kg

kg



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



art. no. 040N
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



art. no. 046
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



art. no. B1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



art. no. 045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering

Additional motor control options



art. no. 090E
Hand switch (low voltage), with magnet on the back to fasten



art. no. 091E
Foot switch (low voltage)



art. no. 092E
1 foot switch fixation mounted at underframe



art. no. 094E
Foot switch rail, lateral, on both sides, security adjustment:
Press the switch rail down = the table is lifted
Lift the switch rail up = the table is lowered



art. no. 096E
Foot switch rail, surround, security adjustment:
Press the switch rail down = the table is lifted
Lift the switch rail up = the table is lowered

Additional safety options



art. no. B1997E
Foot switch cover, made of brushed VA-steel

Paper roll holder



art. no. 010K/015K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



art. no. 010K/F / 015K/F
Paper roll holder, mounted at foot part (not protruding), for width 700 mm / 800 mm

Wall protection bumpers



art. no. B1175 / B1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Head part options



art. no. 1615-EF
Loss of nose opening and upholstery part to close



art. no. 056
Head part 3-sectioned, arm rests downwards adjustable and removable (art. no. 010K/015K not possible)

Side guards



art. no. 1008
2 side guards foldable, RAL 9006 powder coated, made of steel tube, without any breakable plastic components, length: 820 mm

Upholstery options

art. no. 007: Width 800 mm

art. no. B1770 / B1788:
Soft covering for pleasant resting comfort, for width 700 mm/ 800 mm, thickness 80 mm

art. no. B1772 / B1782:
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Serie DX1



3-sectioned versatile treatment tables with hydraulic height adjustment

...for examination, treatment, therapy...

Model 1050-00/H white powder coated (RAL 9010)

Model 1050-03/H white aluminium powder coated (RAL 9006)

Model 1050-04/H silver aluminium powder coated (RAL 9007)

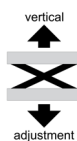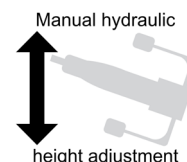

Model 1050-03/H with options:
Movable on twin castors with central lock (046)



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+40°) and downwards (-30°) adjustable by gas spring
- Length head part: 550 mm, integrated nose opening and upholstery part to close
- Fixed middle par
- Length middle part: 480 mm
- Foot part stepless adjustable by gas spring (0° to +50°)
- Length foot part: 920 mm
- Upholstery thickness: 65 mm



Model 1050-03/H with options:
Movable on twin castors with central lock (046)



Technical data:

- Manual foot operated height adjustment: 490 - 940 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width overall: 700 mm
- Length overall: 1950 mm
- Load capacity: 225 kg, according to IEC 60601 tested with fourfold static load
- Load capacity head part: 90 kg


kg



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

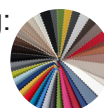

Available as an option (see pricelist for additional options)

Mobility



art. no. 040N
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm




art. no. 046
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



art. no. B1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



art. no. 045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering

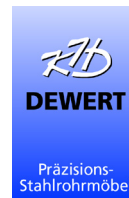

Paper roll holder



art. no. 010K/015K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



art. no. 010K/F / 015K/F
Paper roll holder, mounted at foot part (not protruding), for width 700 mm / 800 mm



art. no. 1008
2 side guards foldable, RAL 9006 powder coated, made of steel tube, without any breakable plastic components, length: 820 mm

Wall protection bumpers




art. no. B1175 / B1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Head part options



art. no. 056
Head part 3-sectioned, arm rests downwards adjustable from -25° to +45° and removable (art. no. 010K/015K not possible)
measures:
Head part: L. 550 mm x W. 380 mm
Arm rests: L. 500 mm x W. 160 mm

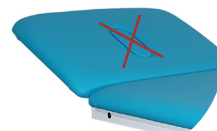

art. no. 1615-EF
Loss of nose opening and upholstery part to close

Upholstery options

art. no. 007
Width 800 mm

art. no. B1770 / B1788
Soft covering for pleasant resting comfort, for width 700 mm/ 800 mm, thickness 80 mm

art. no. B1772 / B1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm



Series DX1


2-sectioned recovery and patient transport table

with hydraulic height adjustment

Model 1090-00/H white powder coated (RAL 9010)

Model 1090-03/H white aluminium powder coated (RAL 9006)

Model 1090-04/H silver aluminium powder coated (RAL 9007)

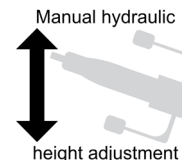

Model 1090-04/H
with option: Storage basket
(B1850)



Basic equipment:

- Pedals for foot-controlled height positioning on both sides of the table
- Movable on twin castors, Ø 100 mm, with central lock, 1 castor with direction lock
- Push handles at head part and foot part, RAL 9006 powder coated
- 2 side guards, foldable, RAL 9006 powder coated, length: 820 mm
- 4 wall protection bumpers, Ø 200 mm
- Upholstery medium-soft (high-quality cold foam cover)
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless adjustable by gas spring (0° - 65°)
- Length head part: 550 mm
- Fixed foot part
- Length foot part: 1400 mm
- Upholstery thickness: 60 mm



Model 1090-04/H
with option: Storage basket (B1850)



Technical data:

- Manual foot operated height adjustment: 500 mm bis 950 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: 700 mm, width overall: 780 mm
- Length of lying surface: 1950 mm, length overall: 2060 mm
- Load capacity: 225 kg, according to IEC 60601 tested with fourfold static load
- Load capacity head part: 90 kg



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



art. no. 1181
Without basic mobility



art. no. B1185
Movable on comfort twin
castors, Ø 125 mm with central lock,
basic height + 30 mm



art. no. B1199
Footbar for central lock,
at foot end

Paper roll holder



art. no. 011K / 013K
Paper roll holder, mounted under head
part, for width 700 mm / 800 mm

Upholstery options

art. no. 016
Width 800 mm

art. no. B1770 / B1780
Soft covering for pleasant resting
comfort, for width 700 mm/ 800 mm,
thickness 80 mm

Add-on elements



art. no. B1410
1 piece of side rail, length: 350 mm



art. no. B1430
1 attachment clamp for side rail,
suitable for pipe: Ø 18 mm



art. no. B1450
Infusion pole, without attachment clamp



art. no. B1850
Storage basket,
measures:
L. 530 x W. 370 x H. 80 mm

NEW


Series DX1

2-sectioned recovery and patient transport table with hydraulic height adjustment

Model 1093-00/H white powder coated (RAL 9010)

Model 1093-03/H white aluminium powder coated (RAL 9006)

Model 1093-04/H silver aluminium powder coated (RAL 9007)

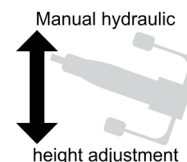

Model 1093-04/H
with basic equipment

Basic equipment:


- Pedals for foot-controlled height positioning on both sides of the table
- Movable on twin castors, Ø 100 mm, with central lock, 1 castor with direction lock
- Push handles at head part and foot part, RAL 9006 powder coated
- 2 side guards, foldable, RAL 9006 powder coated, length: 820 mm
- 4 wall protection bumpers, Ø 200 mm
- Upholstery medium-soft (high-quality cold foam cover)
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless adjustable by gas spring (0° - +75°) (sitting position)
- Length head part: 750 mm
- Fixed foot part
- Length foot part: 1200 mm
- Upholstery thickness: 70 mm

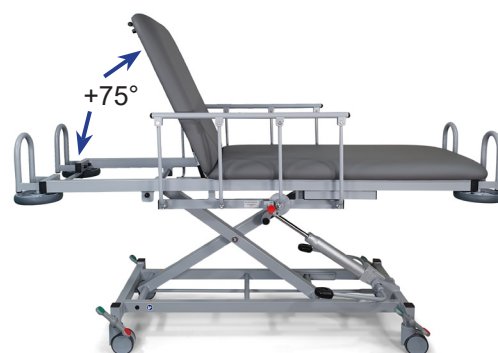

Model 1093-04/H with basic equipment



Technical data:


- Manual foot operated height adjustment: 540 - 990 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: 700 mm, width overall: 780 mm
- Length of lying surface: 1950 mm, length overall: 2060 mm
- Load capacity: 225 kg, according to IEC 60601 tested with fourfold static load
- Load capacity head part: 90 kg



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

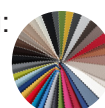

Available as an option (see pricelist for additional options)

Mobility



art. no. 1181
Without basic mobility



art. no. B1185
Movable on comfort twin
castors, Ø 125 mm with central lock,
basic height + 30 mm



art. no. B1199
Footbar for central lock,
at foot end

Paper roll holder



art. no. 011K / 013K
Paper roll holder, mounted under head
part, for width 700 mm / 800 mm

Upholstery options

art. no. 016
Width 800 mm

art. no. B1770 / B1780
Soft covering for pleasant resting
comfort, for width 700 mm/ 800 mm,
thickness 80 mm

Add-on elements



art. no. B1410
1 piece of side rail, length: 350 mm



art. no. B1430
1 attachment clamp for side rail,
suitable for pipe: Ø 18 mm



art. no. B1450
Infusion pole, without attachment clamp



art. no. B1850
Storage basket,
measures:
L. 530 x W. 370 x H. 80 mm




Serie DX1

3-sectioned recovery and patient transport table with hydraulic height adjustment

Model 1095-00/H white powder coated (RAL 9010)

Model 1095-03/H white aluminium powder coated (RAL 9006)

Model 1095-04/H silver aluminium powder coated (RAL 9007)

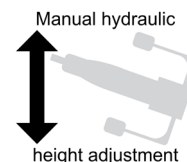

Model 1095-00/H
with options:
Movable on comfort twin castors with central lock (B1185), side rail (B1410)



Basic equipment:

- Pedals for foot-controlled height positioning on both sides of the table
- Movable on twin castors, Ø 100 mm, with central lock, 1 castor with direction lock
- Push handles at head part and foot part, RAL 9006 powder coated
- 2 side guards, foldable, RAL 9006 powder coated, length: 820 mm
- 4 wall protection bumpers, Ø 200 mm
- Upholstery medium-soft (high-quality cold foam cover)
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless adjustable by gas spring (0° - +65°)
- Length head part: 550 mm
- Fixed middle part
- Length middle part: 480 mm
- Foot part stepless adjustable by gas spring (0° to +50°)
- Length foot part: 920 mm
- Upholstery thickness: 70 mm



Model 1095-00/H with options:
Movable on comfort twin castors with central lock (B1185), side rail (B1410)



Technical data:


- Manual foot operated height adjustment: 520 mm bis 970 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: 700 mm, width overall: 780 mm
- Length of lying surface: 1950 mm, length overall: 2060 mm
- Load capacity: 225 kg, according to IEC 60601 tested with fourfold static load
- Load capacity head part: 90 kg



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

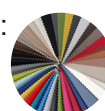

Mobility



art. no. 1181
Without basic mobility



art. no. B1185
Movable on comfort twin
castors, Ø 125 mm with central lock,
basic height + 30 mm



art. no. B1199
Footbar for central lock,
at foot end

Paper roll holder



art. no. 011K / 013K
Paper roll holder, mounted under head
part, for width 700 mm / 800 mm

Upholstery options

art. no. 017
Width 800 mm

art. no. B1770 / B1780
Soft covering for pleasant resting
comfort, for width 700 mm/ 800 mm,
thickness 80 mm

Add-on elements



art. no. B1410
1 piece of side rail, length: 350 mm



art. no. B1430
1 attachment clamp for side rail,
suitable for pipe: Ø 18 mm



art. no. B1450
Infusion pole, without attachment clamp

+ Additional equipment

Cushions



Art.-Nr. 1211-..*
Wedge shaped bolster,
measures (L. x W. x H.):
400 x 400 x 100 - 10 mm

*Add colour number



art. no. 1113-..*
Half roll bolster,
measures (L. x W. x H.):
570 x 200 x 95 mm



art. no. 1111-..*
Neck roll,
measures:
L. 570 x Ø 150 mm



Art.-Nr. 1112-..*
Knee roll,
measures:
L. 570 x Ø 200 mm




Art.-Nr. 1314-..*
Cushion, velcro fixation,
ticking with
polyesterf filling,
measures (L. x W.):
400 x 400 mm

Upholstery protection mat



art. no. 1890
Protection mat, made of PVC,
thickness: 0.8 mm,
measures: 700 x 600 mm

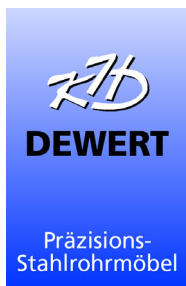

K.H. Dewert GmbH
Vollmestraße 7
33649 Bielefeld
Germany

Post office 14 40 28
33634 Bielefeld
Germany

Phone + 49 (0) 5 21 / 4 00 27-0
Fax + 49 (0) 5 21 / 4 00 27-27

Internet: www.khdewert.de
E-mail: info@khdewert.de
YouTube: K.H. Dewert

You will find the whole Dewert product range with additional articles and equipment on our website and in our further catalogues. We look forward to hear from you!



EXAMINATION AND THERAPY TABLES

SERIES DX1

2021.2_3

