

# Simple and clean design

**Day surgery table S** is designed to give you **maximum support** during day surgery procedures. It allows the surgeon and his team to **access** the patient as **easily** as possible. The table top is **65 cm wide** and features a narrow head section that includes a a breathing hole and a plug as part of standard equipment. The **padding** is **8 cm thick** and upholstered in **antistatic black imitation leather**. It is **sealed** on the underside with welded and taped seams, which prevent the ingress of liquids into the interior. This makes cleaning and ensuring hygiene standards even easier.

## Reliable

The **four-section table** with the enviable **safe working load of 250 kg** and sturdy construction is an ultimately reliable piece of equipment. The built-in software prevents the table from hitting the floor when adjusting its height. In the lowest position, the table stops just above the floor level. The Day surgery table comes with an **additional battery backup** so that in case of a power outage you can reset the position of the table and complete your tasks.

## Flexible

The electric table features a **circular switch** for height adjustment. This is crucial during procedures, as the surgeon can **adjust the height** regardless of where they stand **hands-free**. The middle two sections of the table are electrically supported which makes it easier to position your patients. The **hand switch** is simple to use and features an **innovative memory function**. You can save two positions of your choice in addition to a completely flat position. The head and leg sections are supported by gas springs and adjust manually. The table is standardly equipped with **fixation rails** and **armrests** that are easily adjustable.



[in](#) [f](#) Bring treatment to a higher level!

Novak M d.o.o.

Pod smrekami 5  
SI - 1218 Komenda

E: [info@novak-m.si](mailto:info@novak-m.si)  
T: +386 1 5620117


[www.novak-m.com](http://www.novak-m.com)



# Day surgery




TOP VIEW

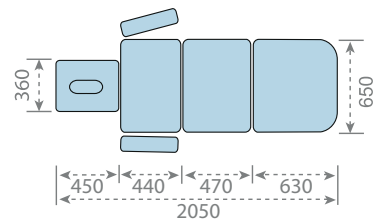

# DAY SURGERY TABLE S

Model: **ETSOP-65**

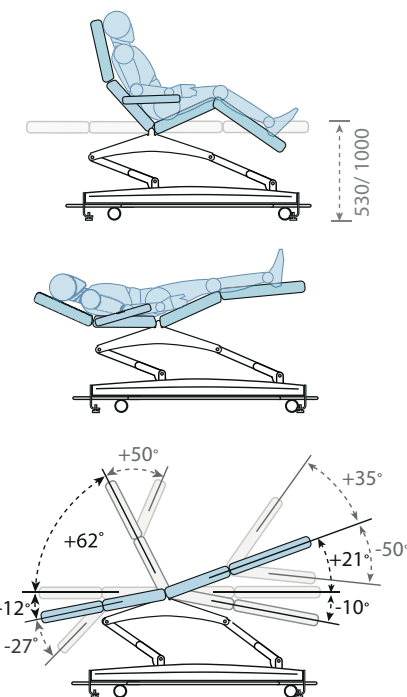

POSITIONS

## Narrow head section

- with a breathing hole and a plug
- manually adjustable

## Adjustable armrests

- attached to fixation rails
- adjustable height and incline

## Black antistatic upholstery

- 65 cm wide
- seamless and fully sealed
- 8 cm thick padding with a memory foam layer
- wider lying surface due to the padding finished at a right angle

## Table base

- ABS plastic cover is screwless for easier cleaning
- increased stability due to standard 250 kg safe working load

## Fixation rails

- along back and sitting sections

## Soft retractable system

- antistatic castors
- includes a shock absorber to soften the release of castors

## Circular switch

- simple up & down operation from all sides of the table

## Battery back-up

- built-in for case of power failure



## Hand switch

- for adjustment of back and sitting section
- with innovative memory function

Equipment and accessories	ETSOP-65
safe working load 250 kg	●
circular switch for electric height adjustment	●
antistatic, seamless and sealed upholstery with memory foam (8 cm)	●
hand switch for electric adjustment of height, back and sitting section	●
built-in battery backup in case of power failure	●
adjustable armrests on fixation rails	●
fixation rails at back and sitting section	●
antistatic retractable castors with soft release (Ø 75 mm)	●
plug connector for potential equalization – PPEQ	○
safety switch – EVV-S	○

standard ● optional ○